




ANALYTICAL REPORT

4-Chloropentedrone (C₁₂H₁₆ClNO)

1-(4-chlorophenyl)-2-(methylamino)pentan-1-one

Remark – other active cpd. detected: **none**

Sample ID:	1708-16
Sample description:	powder - white
Sample type:	RM-reference material
Comments ¹ :	CAY Lot#04829156; RESPONSE -purchasing
Date of entry:	1/6/2017

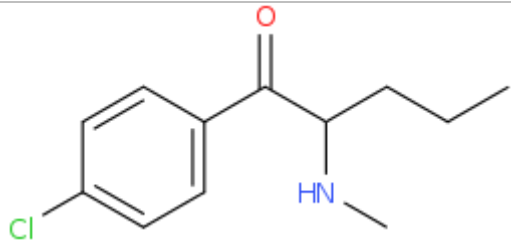Substance identified-structure ² (base form)	
Systematic name:	1-(4-chlorophenyl)-2-(methylamino)pentan-1-one
Other names:	4-CPD; 4-chloro- α -methylamino-Pentiophenone; 4-chloro- α -methylamino-Valerophenone
Formula (per base form)	C ₁₂ H ₁₆ ClNO
M _w (g/mol)	225,72
Salt form:	HCl
StdInChIKey (for base form)	OPCMVVKRCLOEDQ-UHFFFAOYSA-N
Other active cpd. detected	none
Add.info (purity..)	98%

¹ This report has been produced with the financial support of the Prevention of and fight against crime Programme of the European Union (grant agreement number JUST/2013/ISEC/DRUGS/AG/6413). The contents of this report are the sole responsibility of the National Forensic Laboratory and can in no way be taken to reflect the views of the European Commission.

² Created by OPSIN free tool: <http://opsin.ch.cam.ac.uk/> DOI: 10.1021/ci100384d



Report updates

date	comments (explanation)

Supporting information

Analytical technique:	applied	remarks
GC-MS (EI ionization)	+	NFL GC-RT (min): 4,86 BP(1): 86; BP(2): 44,BP(3) :111,
FTIR-ATR	+	direct measurement
GC-IR (condensed phase)	+	always as base form

1. GC-MS (Agilent): GC-method is RT locked to tetracosane (9.258 min). Injection volume 1 ml and split mode (1:50). Injector temperature: 280 °C. Chromatographic separation: on column HP1-MS (100% dimethylpolysiloxane), length 30 m, internal diameter 0.25 mm, film thickness 0.25 µm. Carrier gas He: flow-rate 1.2 ml/min. GC oven program: 170 °C for 1 min, followed by heating up to 190 °C at rate 8 °C/min, then heating up to 293 °C at a rate of 18 °C/min, hold for 7.1 min, then heating at 50 °C/min up to 325 °C and finally 6.1 min isothermal. MSD source EI = 70 eV. GC-MS transfer line T= 235°C, source and quadropole temperatures 280°C and 180°C, respectively. Scan range m/z scan range: from 50 (30 until 6 min.) to 550 (300 until 6 min) amu.

2. FTIR-ATR (Perkin Elmer): scan range 4000-400 cm⁻¹; resolution 4cm⁻¹

3. GC- (MS)-IR condensed phase (GC-MS (Agilent) & IR (Spectra analyses-Danny)

GC-method: Injection volume 1 ml and split mode (1:5). Injector temperature 280 °C. Chromatographic separation as above (1). Split MS : IR = 1 : 9.

MSD source EI = 70 eV. GC-MS transfer line T= 235°C, source and quadropole temperatures 280°C and 180°C, respectively. Scan range m/z scan range: from 50 (30 until 6 min.) to 550 (300) amu.


IR (condensed (solid) phase): IR scan range 4000 to 650, resolution 4 cm⁻¹.

4. HPLC-TOF for exact monoisotopic mass and empirical formula control - results are not shown in the report.


FIGURES OF SPECTRA

MS (EI)

Abundance

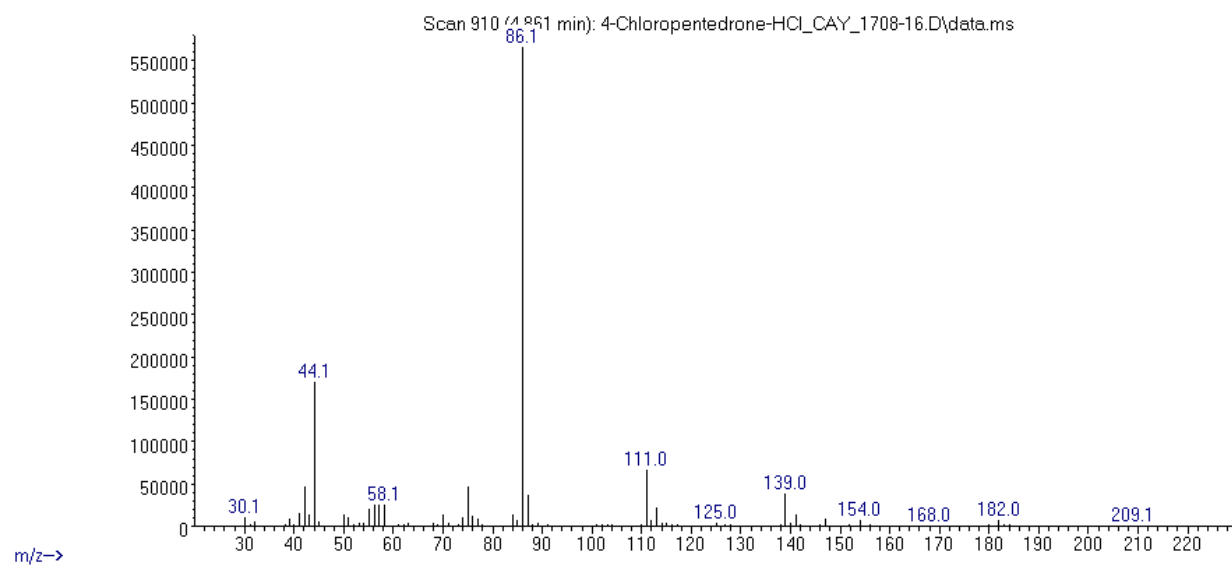


FTIR-ATR (sample as received)

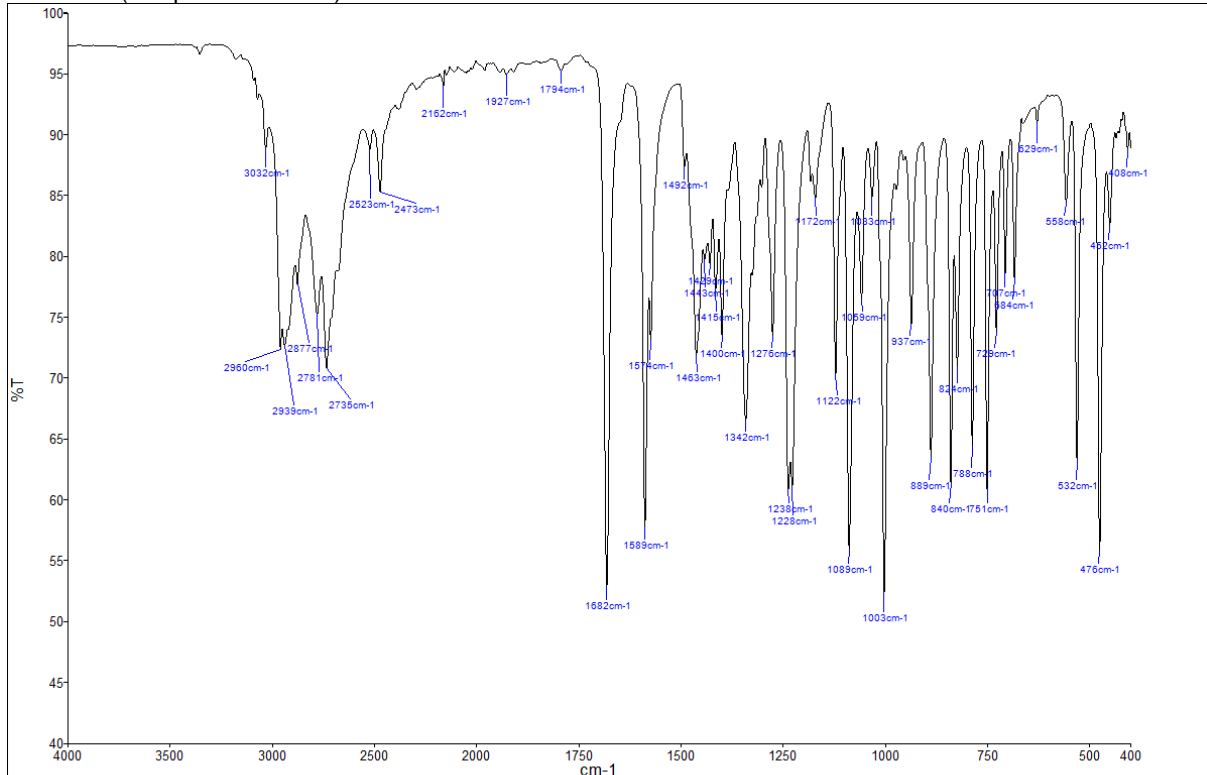

Name	Description
4-Chloropentredone-HCl_1708-16_CAY_IR-A	ID:1708-16;1-(4-chlorophenyl)-2-(methylamino)pentan-1-one-HCl_Lot#0482915-6

IR-Condensed phase

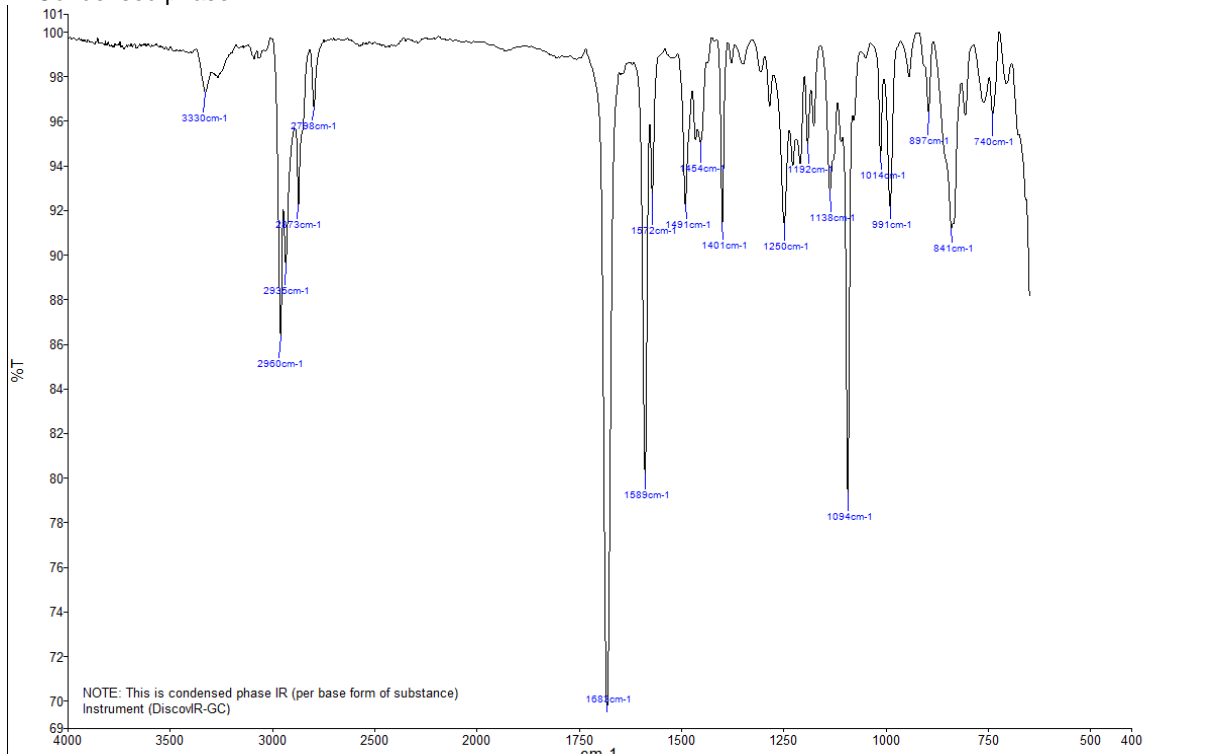

NOTE: This is condensed phase IR (per base form of substance)
Instrument (DiscovIR-GC)

Name	Description
4-Chloropentredone-HCl_1708-16_Cay_IR-C.spc	Sample_ID:1708-16;1-(4-chlorophenyl)-2-(methylamino)pentan-1-one-HCl_Lot#0482915-6