



NACIONALNI FORENZIČNI LABORATORIJ
NATIONAL FORENSIC LABORATORY

Vodovodna 95
1000 Ljubljana
SLOVENIJA

T: +386 (0)1 428 44 93

E: nfl@policija.si

www.policija.si

ANALYTICAL REPORT

2C-P-NBOMe (C₂₁H₂₉N₃O₃)

[2-(2,5-dimethoxy-4-propylphenyl)ethyl][(2-methoxyphenyl)methyl]amine

Remark – other active cpd. detected: **none**

Sample ID:	1681-16
Sample description:	crystallinic - white
Sample type:	RM-reference material
Comments ¹ :	CAY Lot#0459579; I-SEE project RM and report
Date of entry:	10/13/2016

Substance identified- structure ² (base form)	
Systematic name:	[2-(2,5-dimethoxy-4-propylphenyl)ethyl][(2-methoxyphenyl)methyl]amine
Other names:	25P-NB2OMe; 2-(4-Propyl-2,5-dimethoxyphenyl)-N-(2-methoxybenzyl)ethanamine
Formula (per base form)	C ₂₁ H ₂₉ N ₃ O ₃
M _w (g/mol)	344,47
Salt form:	HCl
StdInChIKey	QHEBYIJRKGDGFGH-UHFFFAOYSA-N
Compound Class	Phenethylamines
Other active cpd. detected	none
Add.info (purity..)	>98.5 % (as HCL)

¹ This report has been produced with the financial support of the Prevention of and fight against crime Programme of the European Union (grant agreement number JUST/2013/ISEC/DRUGS/AG/6426). The contents of this report are the sole responsibility of the National Forensic Laboratory and can in no way be taken to reflect the views of the European Commission.

² Created by OPSIN free tool: <http://opsin.ch.cam.ac.uk/> DOI: 10.1021/ci100384d

Report updates

date	comments (explanation)

Supporting information

Analytical technique:	applied	remarks
GC-MS (EI ionization)	+	NFL GC-RT (min): 10,09 BP(1): 121; BP(2): 150,BP(3) :91,
FTIR-ATR	+	
GC-IR (condensed phase)	+	

1. GC-MS (Agilent): GC-method is RT locked to tetracosane (9.258 min). Injection volume 1 ml and split mode (1:50). Injector temperature: 280 °C. Chromatographic separation: on column HP1-MS (100% dimethylpolysiloxane), length 30 m, internal diameter 0.25 mm, film thickness 0.25 µm. Carrier gas He: flow-rate 1.2 ml/min. GC oven program: 170 °C for 1 min, followed by heating up to 190 °C at rate 8 °C/min, then heating up to 293 °C at a rate of 18 °C/min, hold for 6.1 min, then heating at 50 °C/min up to 325 °C and finally 6.1 min isothermal. MSD source EI = 70 eV. GC-MS transfer line T= 235°C, source and quadropole temperatures 280°C and 180°C, respectively. Scan range m/z scan range: from 50 (30 until 6 min.) to 550 (300 until 6 min) amu.

2. FTIR-ATR (Perkin Elmer): scan range 4000-400 cm⁻¹; resolution 4cm⁻¹

3. GC- (MS)-IR condensed phase (GC-MS (Agilent) & IR (Spectra analyses-Danny)

GC-method: Injection volume 1 ml and split mode (1:5). Injector temperature 280 °C. Chromatographic separation as above **(1)**. Split MS : IR = 1 : 9.


MSD source EI = 70 eV. GC-MS transfer line T= 235°C, source and quadropole temperatures 280°C and 180°C, respectively. Scan range m/z scan range: from 50 (30 until 6 min.) to 550 (300) amu.

IR (condensed (solid) phase): IR scan range 4000 to 650, resolution 4 cm⁻¹.


FIGURES OF SPECTRA

MS (EI)

Abundance

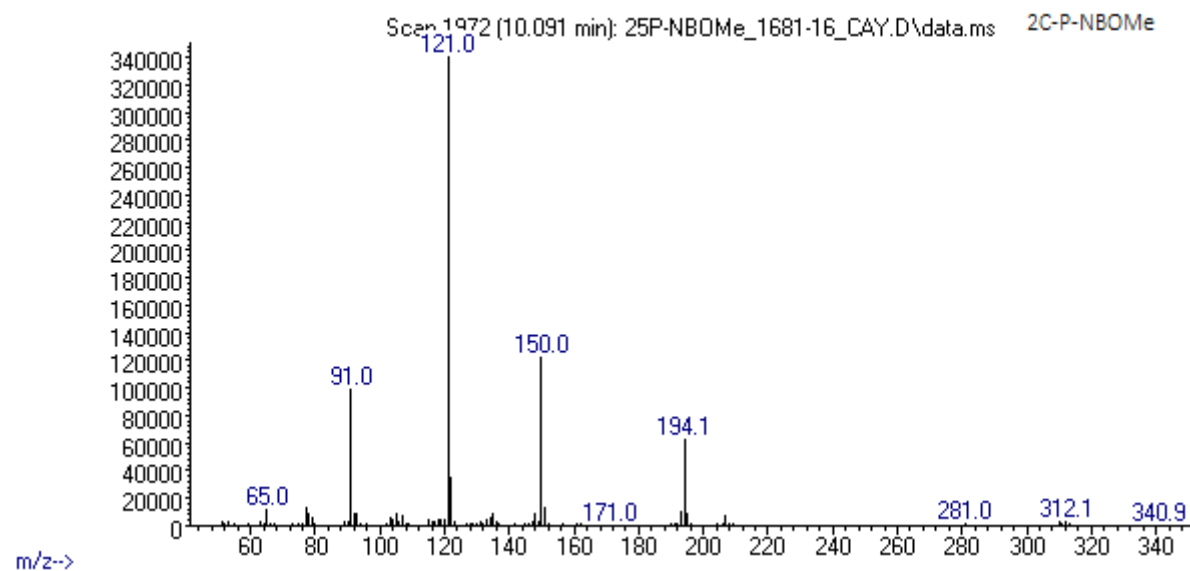

FTIR-ATR (direct measurement)

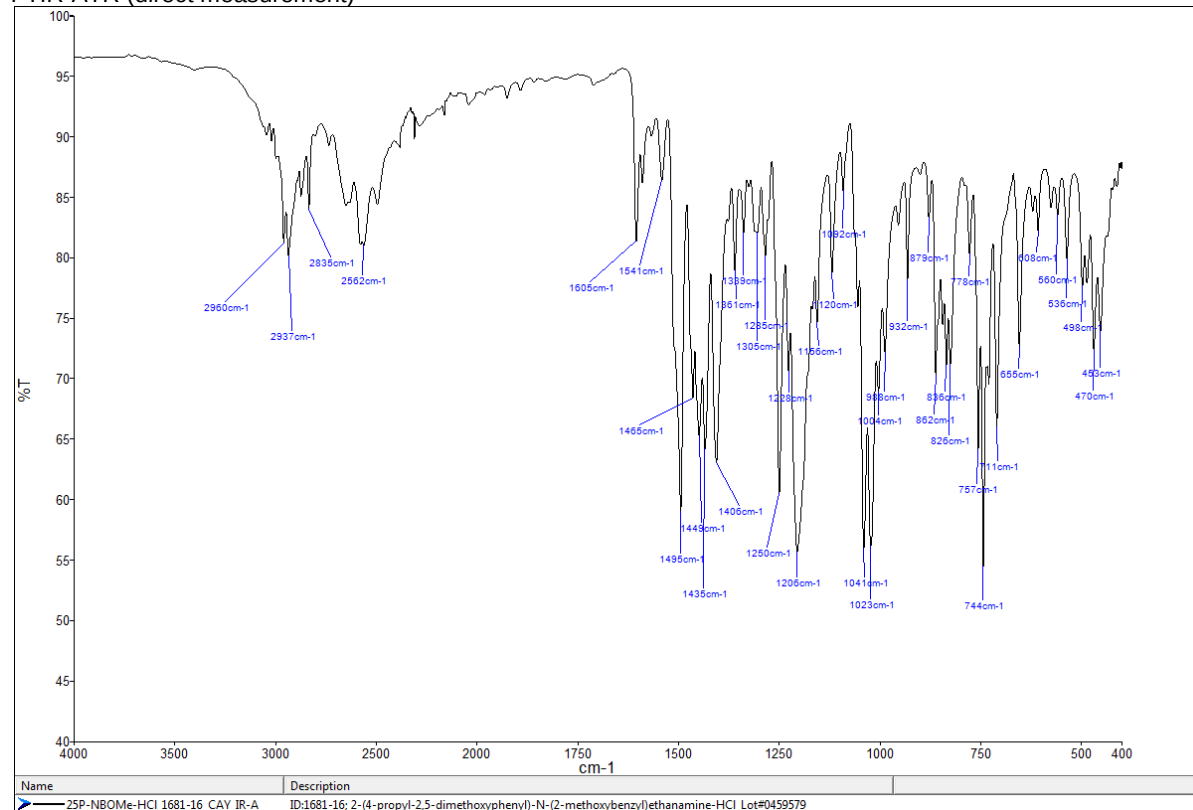

GC-IR-Condensed phase

