




ANALYTICAL REPORT

3-Hydroxyphenazepam (C₁₅H₁₀BrClN₂O₂)

7-Bromo-5-(2-chlorophenyl)-3-hydroxy-1,3-dihydro-2H-1,4-benzodiazepin-2-one

Remark – other active cpd. detected: **none**

Sample ID:	1235-15
Sample description:	powder - light grey
Sample type:	RM-reference material
Comments ¹ :	Chiron AS Lot#14203; RESPONSE -purchasing
Date of entry:	8/31/2015

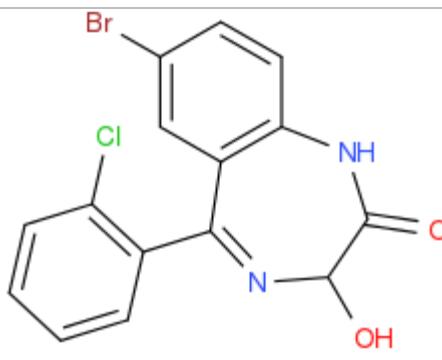Substance identified-structure ² (base form)	
Systematic name:	7-Bromo-5-(2-chlorophenyl)-3-hydroxy-1,3-dihydro-2H-1,4-benzodiazepin-2-one
Other names:	
Formula (per base form)	C ₁₅ H ₁₀ BrClN ₂ O ₂
M _w (g/mol)	365,61
Salt form:	base
StdInChIKey	KRJJKUWAZOWXNV-UHFFFAOYSA-N
Compound Class	Benzodiazepines
Other active cpd. detected	none
Add.info (purity..)	98%

¹ This report has been produced with the financial support of the Prevention of and fight against crime Programme of the European Union (grant agreement number JUST/2013/ISEC/DRUGS/AG/6413). The contents of this report are the sole responsibility of the National Forensic Laboratory and can in no way be taken to reflect the views of the European Commission.

² Created by OPSIN free tool: <http://opsin.ch.cam.ac.uk/> DOI: 10.1021/ci100384d



Report updates

date	comments (explanation)

Supporting information

Analytical technique:	applied	remarks
GC-MS (EI ionization)	+	NFL GC-RT (min): 9,84 BP(1): 320; BP(2): 318,BP(3) :75,
FTIR-ATR	+	direct measurement
GC-IR (condensed phase)	+	spectrum is always for the base form of compound

GC-MS (Agilent):

GC-method is RT locked to tetracosane (RT=9.53 min).

Injection volume 1 ml and split mode (1:50) .

Injector temperature: 280 °C.

Chromatographic separation

Column: HP1-MS (100% dimethylpolysiloxane), length 30 m, internal diameter 0.25 mm, film thickness 0.25 mm.

Carrier gas He: flow-rate 1.2 ml/min. GC oven program: 170 °C for 1 min, followed by heating up to 293 °C at a rate of 18 °C/min, hold for 6.1 min, then heating at 50 °C/min up to 325 °C and finally 2.8 min isothermal.

MSD source EI = 70 eV. GC-MS transfer line T= 235°C, source and quadrupole temperatures 280°C and 180°C, respectively. Scan range m/z scan range: from 50 (40) to 550 amu.

FTIR-ATR (Perkin Elmer): scan range 4000-400 cm⁻¹; resolution 4cm⁻¹

GC- (MS)-IR condensed phase (GC-MS (Agilent) & IR (Spectra analyses-Danny) IR scan range 4000 to 700, resolution 4cm⁻¹

GC-method:

Injection volume 1 ml and split mode (1:5) .

Injector temperature: 280 °C.

Chromatographic separation

Column: HP1-MS (100% dimethylpolysiloxane), length 30 m, internal diameter 0.25 mm, film thickness 0.25 mm.

Carrier gas He: flow-rate 1.2 ml/min. GC oven program: 170 °C for 1 min, followed by heating up to 293 °C at a rate of 18 °C/min, hold for 6.1 min, then heating at 50 °C/min up to 325 °C and finally 2.8 min isothermal.

Split MS : IR : (1:9)


MSD source EI = 70 eV. GC-MS transfer line T= 235°C, source and quadrupole temperatures 280°C and 180°C, respectively. Scan range m/z scan range: from 50 (40) to 550 amu.

IR (condensed phase): IR scan range 4000 to 700, resolution 4cm⁻¹


FIGURES OF SPECTRA

MS (EI)

Abundance

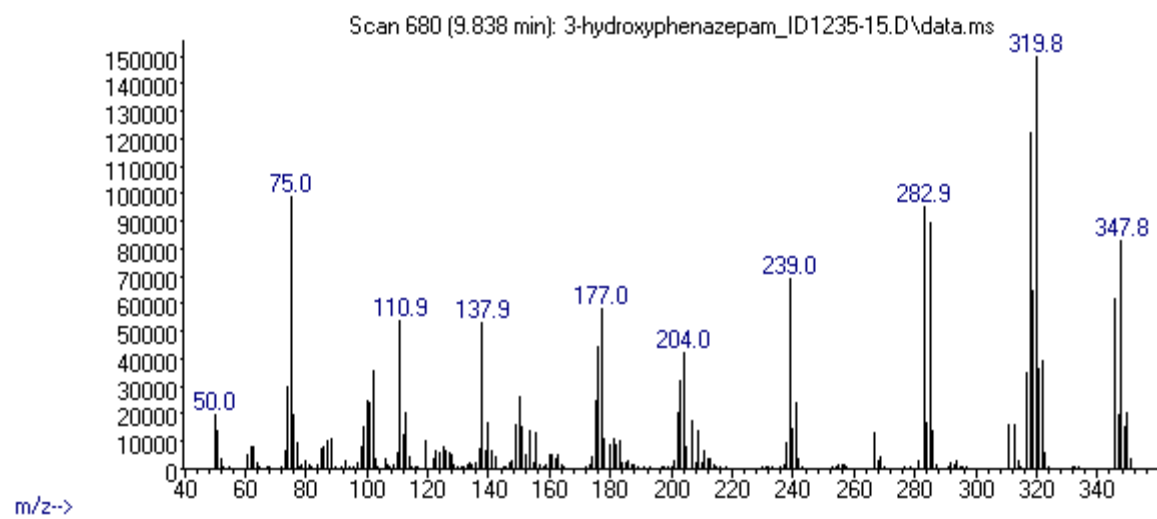

FTIR-ATR

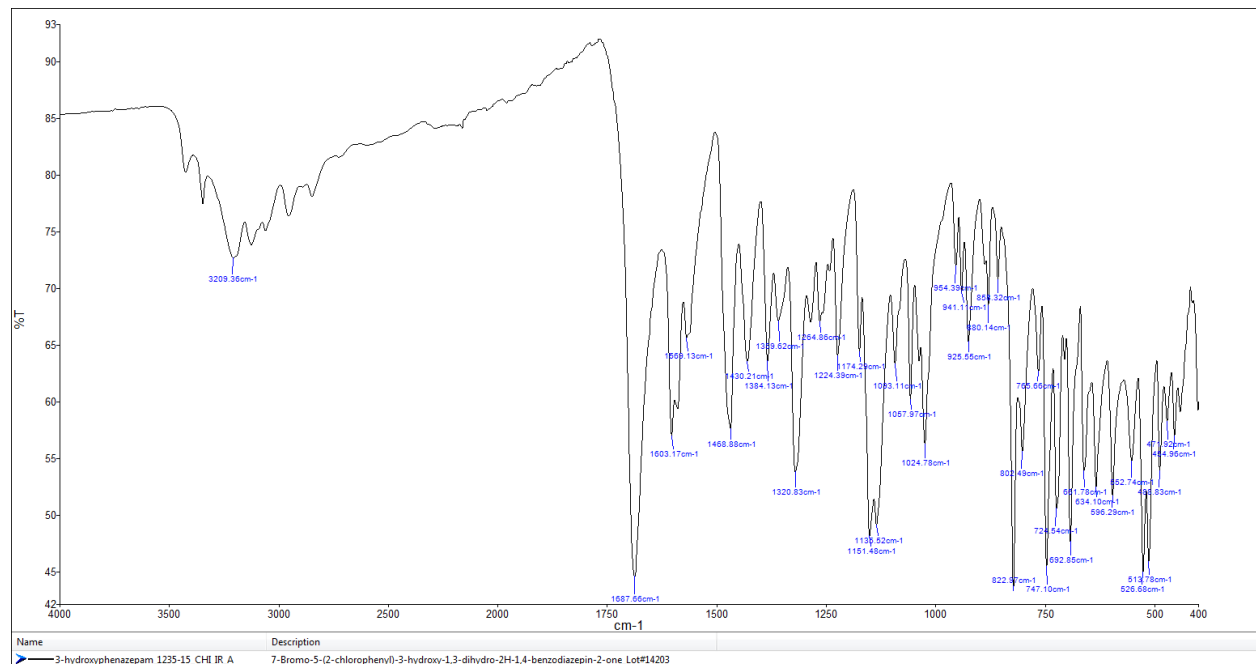

IR-condensed phase

