



NACIONALNI FORENZIČNI LABORATORIJ
NATIONAL FORENSIC LABORATORY

Vodovodna 95
1000 Ljubljana
SLOVENIJA

T: +386 (0)1 428 44 93

E: nfl@policija.si

www.policija.si


ANALYTICAL REPORT

3F-Phenetrazine (C₁₂H₁₆FNO)

3-ethyl-2-(3-fluorophenyl)morpholine

Remark – other NPS detected: **none**

Sample ID:	1292-15
Sample description:	powder - white
Sample type:	test purchase /RESPONSE -purchasing
Comments ¹ :	
Date of entry into NFL database: http://www.policija.si/apps/nfl_response_web/seznam.php	10/14/2015

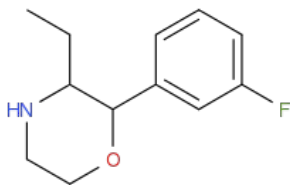Substance identified- structure ² (base form)	
Systematic name	3-ethyl-2-(3-fluorophenyl)morpholine
Other names	
Formula (per base form)	C ₁₂ H ₁₆ FNO
M _w (g/mol)	209.26
Salt form	HCl
StdInChIKey	GPTHUZORYIVVNJ-UHFFFAOYSA-N
Compound Class	Others
Other NPS detected	none
Add.info (purity..)	pure by GC-MS, HPLC-TOF,

¹ This report has been produced with the financial support of the Prevention of and fight against crime Programme of the European Union (grant agreement number JUST/2013/ISEC/DRUGS/AG/6413). The contents of this report are the sole responsibility of the National Forensic Laboratory and can in no way be taken to reflect the views of the European Commission.

² Created by OPSIN free tool: <http://opsin.ch.cam.ac.uk/> DOI: 10.1021/ci100384d

Report updates

date	comments (explanation)
25/10/2016	Figure of IR condensed spectrum replaced.

Instrumental methods (if applied) in NFL

1. GC-MS (Agilent): GC-method is RT locked to tetracosane (RT=9.53 min). Injection volume 1 ml and split mode (1:50). Injector temperature: 280 °C. Chromatographic separation: on column HP1-MS (100% dimethylpolysiloxane), length 30 m, internal diameter 0.25 mm, film thickness 0.25 mm. Carrier gas He: flow-rate 1.2 ml/min. GC oven program: 170 °C for 1 min, followed by heating up to 293 °C at a rate of 18 °C/min, hold for 6.1 min, then heating at 50 °C/min up to 325 °C and finally 2.8 min isothermal. MSD source EI = 70 eV. GC-MS transfer line T= 235 °C, source and quadropole temperatures 280 °C and 180 °C, respectively. Scan range m/z scan range: from 50 (40) to 550 amu.

2. HPLC-TOF (Agilent): 6230B TOF with Agilent 1260 Infinity HPLC with binary pump, column: Zorbax Eclipse XDB-C18, 50 x 4.6 mm, 1.8 micron. Mobile phases (A) 0.1% formic acid and 1mM ammonium formate in water; (B) 0.1% formic acid in methanol (B). Gradient: starting at 5% B, changing to 40% B over 4 min, then to 70% over 2 min and in 5 min to 100%, hold 1 min and back to 5%, equilibration for 1.7 min. The flow rate: 1.0 ml/min; Injection volume 1 µl. MS parameters: 2GHz, Extended Dynamic range mode to a maximum of 1700 amu, acquisition rate 1.30 spectra/sec. Sample ionisation: by Agilent Jet Stream technology (Dual AJS ESI). Ion source: positive ion scan mode with mass scanning from 82 to 1000 amu. Other TOF parameters: drying gas (N2) and sheath temperature 325 °C; drying gas flow rate 6 l/min; sheath gas flow rate 8 l/min; nebulizer 25 psig; Vcap. 4000 V; nozzle 2000 V; skimmer 65 V; fragmentor 175 V and Octopole RF 750 V.

3. FTIR-ATR (Perkin Elmer): scan range 4000-400 cm⁻¹; resolution 4cm⁻¹

4. GC- (MS)-IR condensed phase (GC-MS (Agilent) & IR (Spectra analyses-Danny))
GC-method: Injection volume 1 ml and split mode (1:5). Injector temperature 280 °C. Chromatographic separation as above **(1)**. Split MS : IR = 1:9.
MSD source EI = 70 eV. GC-MS transfer line T= 235 °C, source and quadropole temperatures 280 °C and 180 °C, respectively. Scan range m/z scan range: from 50 (40) to 550 amu.
IR (condensed phase): IR scan range 4000 to 650, resolution 4 cm⁻¹.

Supporting information


Solubility in	result/remark
CH ₂ Cl ₂	soluble
MeOH	soluble
H ₂ O	soluble

Analytical technique:	applied	remarks
GC-MS (EI ionization)	+	NFL GC-RT (min): 3.96 BP(1): 56; BP(2): 85,BP(3) :70,
HPLC-TOF	+	Exact mass (theoretical): 209.1216; measured value Δppm:-1.8; formula:C12H16FNO
FTIR-ATR	+	direct measurement
FTIR (condensed phase) always as base form	+	
NMR		
validation		
other		


ANALYTICAL RESULTS

MS (EI)


Abundance

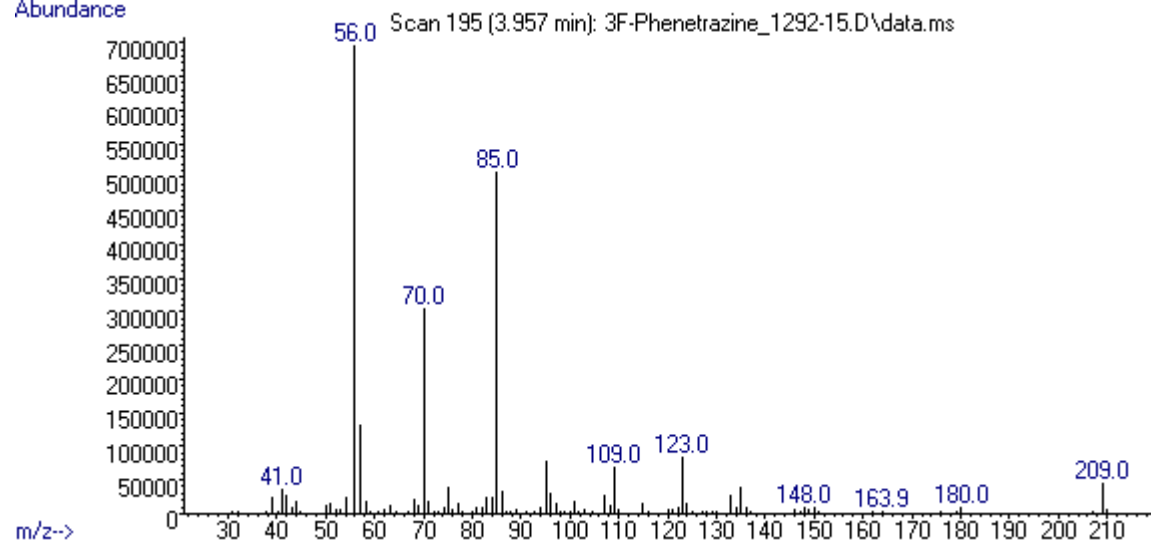

FTIR-ATR - direct measurement

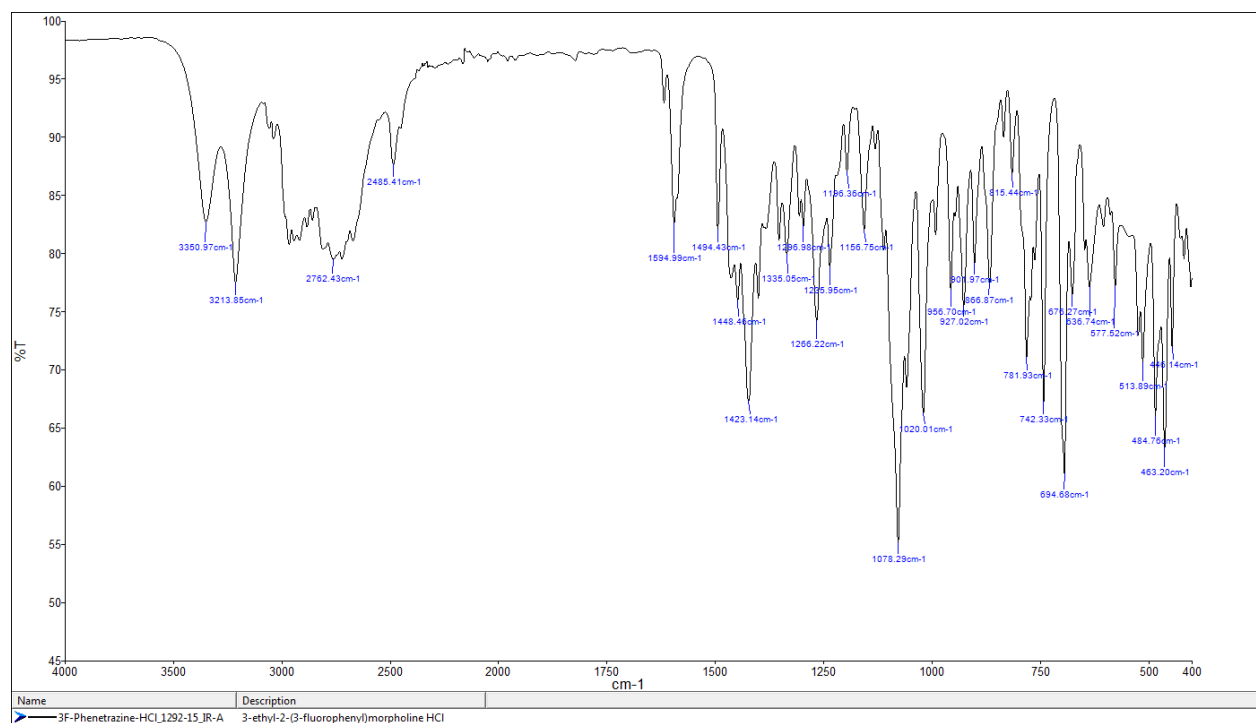


IR (condensed phase)

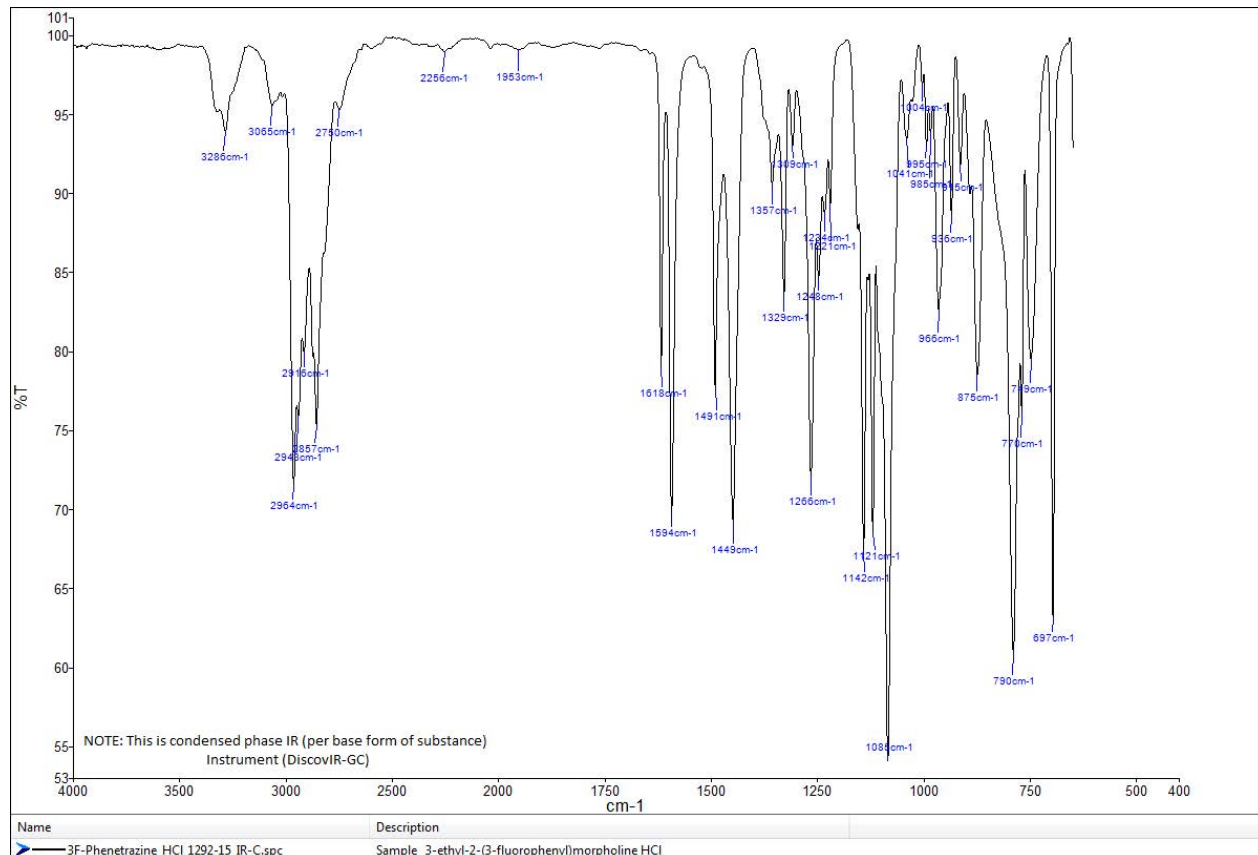

Peak Integration Report

Sample Name:	3F-Phenetrazine_1292-15_IC	Inj. Vol.:	25,00
Injection Type:	Unknown	Dilution Factor:	1,0000
Program:	ANIONI	Operator:	kemija
Inj. Date / Time:	30-sep-2015 / 12:15	Run Time:	41,99

No.	Time min	Peak Name	Peak Type	Area $\mu\text{S}\cdot\text{min}$	Height μS	Amount mg/L
1,00	9,71	Chloride	BMB	4,20	17,25	n.a.
		TOTAL:		4,20	17,25	0,00

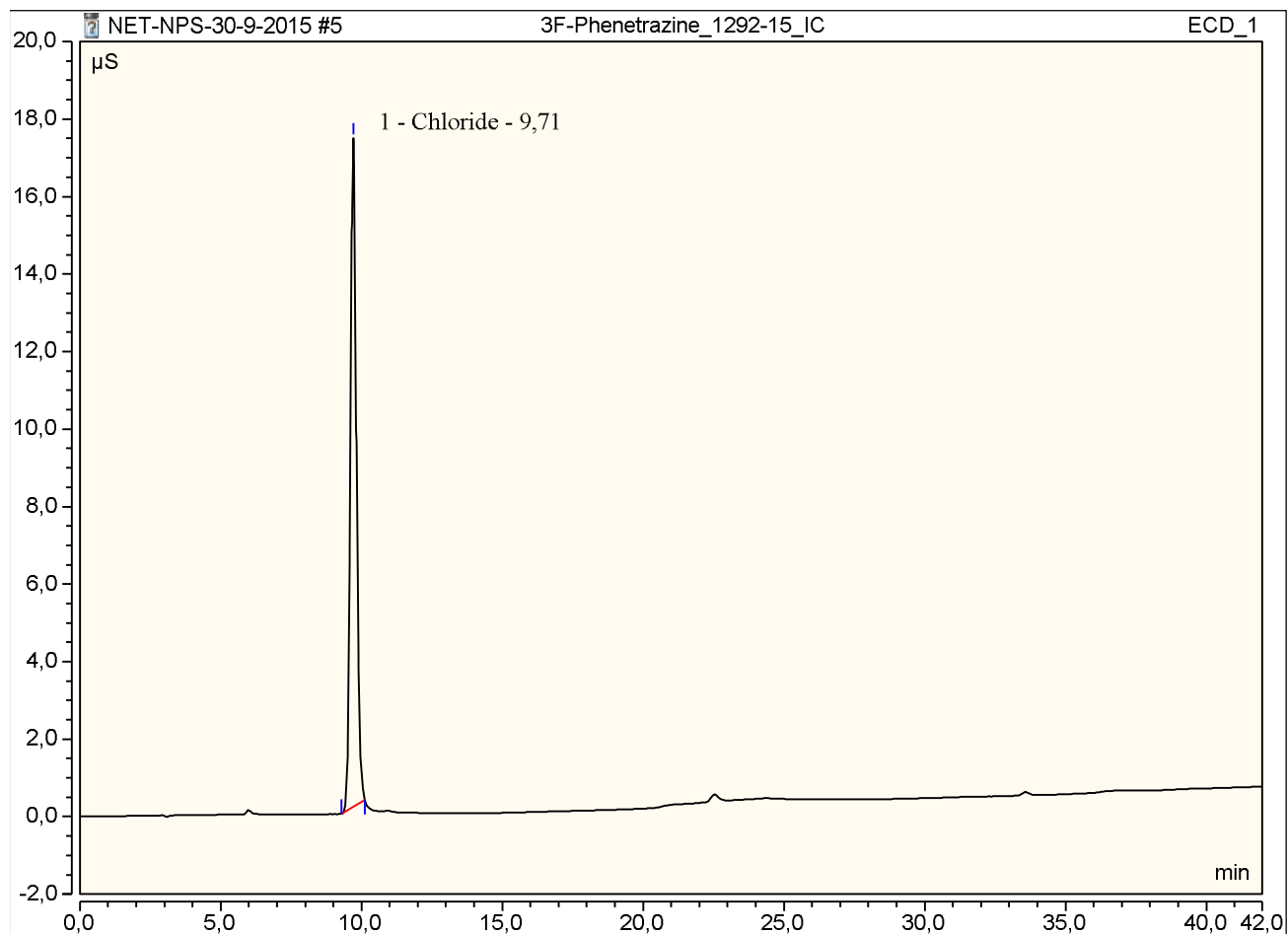

TOF REPORT


Data File	3F-Phenetrazine_1292-15_TOF.d	Sample Name	3F-Phenetrazine
Sample Type	Sample	Position	P2-B7
Instrument Name	6230B TOF LC-MS	User Name	TG
Acq Method	general-28052015-XDB-C18-ESI-poz.m	Acquired Time	9/28/2015 1:20:52 PM
IRM Calibration Status	Success	DA Method	Droge_Default.m
Comment	extract in MeOH		

Compound Table


Label	Compound Name	Obs. RT	Obs. Mass
Cpd 1: 3F-Phenetrazine	3F-Phenetrazine	4.858	209.122

Name	Obs. m/z	Obs. RT	Obs. Mass	DB RT	DB Formula	DB Mass	DB Mass Error (ppm)
3F-Phenetrazine	210.1292	4.858	209.122	4.86	C12 H16 F N O	209.1216	-1.8


Compound Chromatograms

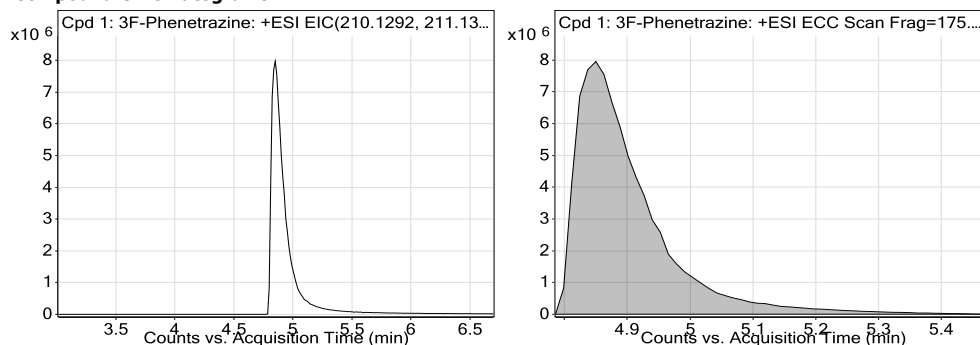

MFE MS Zoomed Spectrum

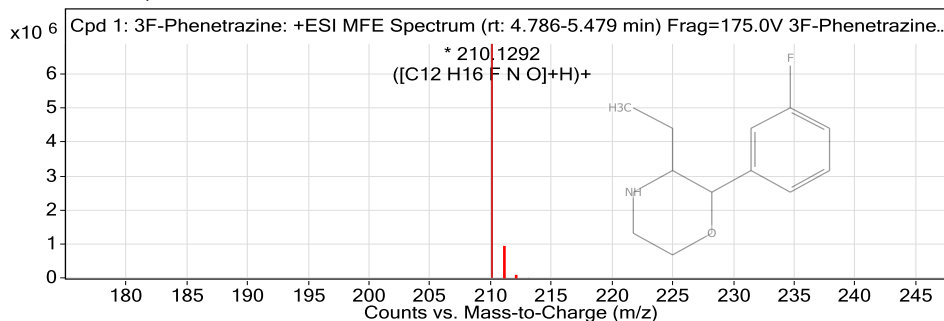

MS Zoomed Spectrum

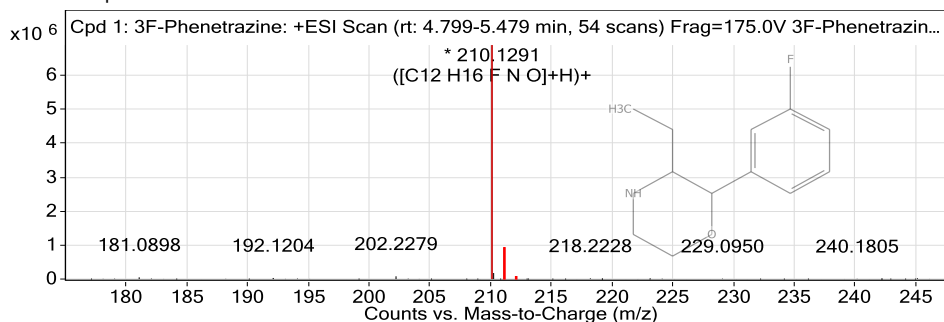


MS Spectrum Peak List

Obs. m/z	Charge	Abund	Formula	Ion/Isotope
210.1292	1	6889866	C12 H16 F N O	(M+H)+
211.133	1	904715.8	C12 H16 F N O	(M+H)+
212.1355	1	64936.09	C12 H16 F N O	(M+H)+
213.1385	1	3527.29	C12 H16 F N O	(M+H)+


--- End Of Report ---



REPORT

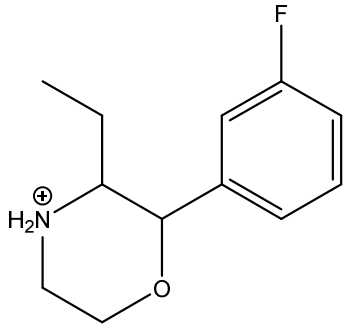Sample ID:	1292-15
Our notebook code:	P-1292-15
NMR sample preparation:	15 mg dissolved in 0.7 mL DMSO- d_6
NMR experiments:	^1H , ^{13}C .
Proposed structure:	
Chemical name:	3-ethyl-2-(3-fluorophenyl)morpholin-4-ium
Comments:	- Structure elucidation based on 1D NMR spectra - Compound is pure by NMR.
Supporting information:	Copies of ^1H and ^{13}C NMR spectra
Author:	Prof. Dr. Janez Košmrlj, Doc. Dr. Krištof Kranjc
Date of report:	December 8, 2015

P-1292-15

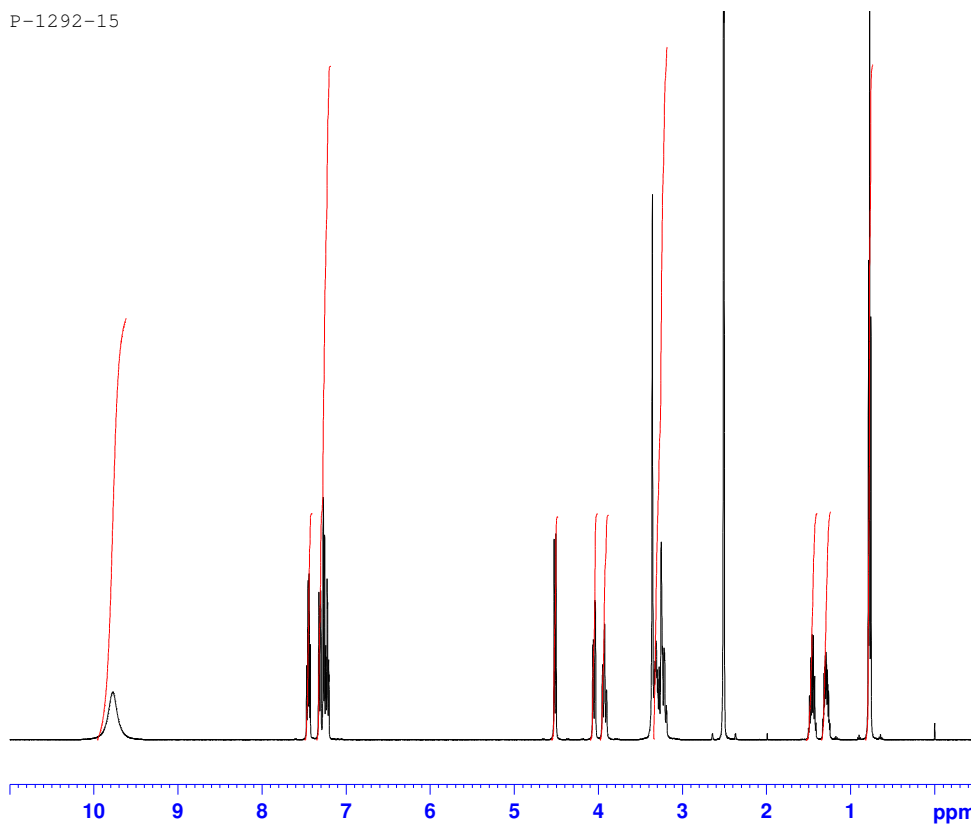


Current Data Parameters
NAME P-1292-15
EXPNO 1
PROCNO 1

F2 - Acquisition Parameters
Date_ 20151128
Time 21.51
INSTRUM spect
PROBHD 5 mm PABBO BB-
PULPROG zg30
TD 65536
SOLVENT DMSO
NS 16
DS 2
SWH 10000.000 Hz
FIDRES 0.152588 Hz
AQ 3.2768500 sec
RG 90.5
DW 50.000 usec
DE 6.50 usec
TE 296.0 K
D1 1.00000000 sec
TD0 1

===== CHANNEL f1 =====
SFO1 500.1330885 MHz
NUC1 1H
P1 8.90 usec
PLW1 26.00000000 W

F2 - Processing parameters
SI 65536
SF 500.1300002 MHz
WDW EM
SSB 0
LB 0.30 Hz
GB 0
PC 1.00

P-1292-15

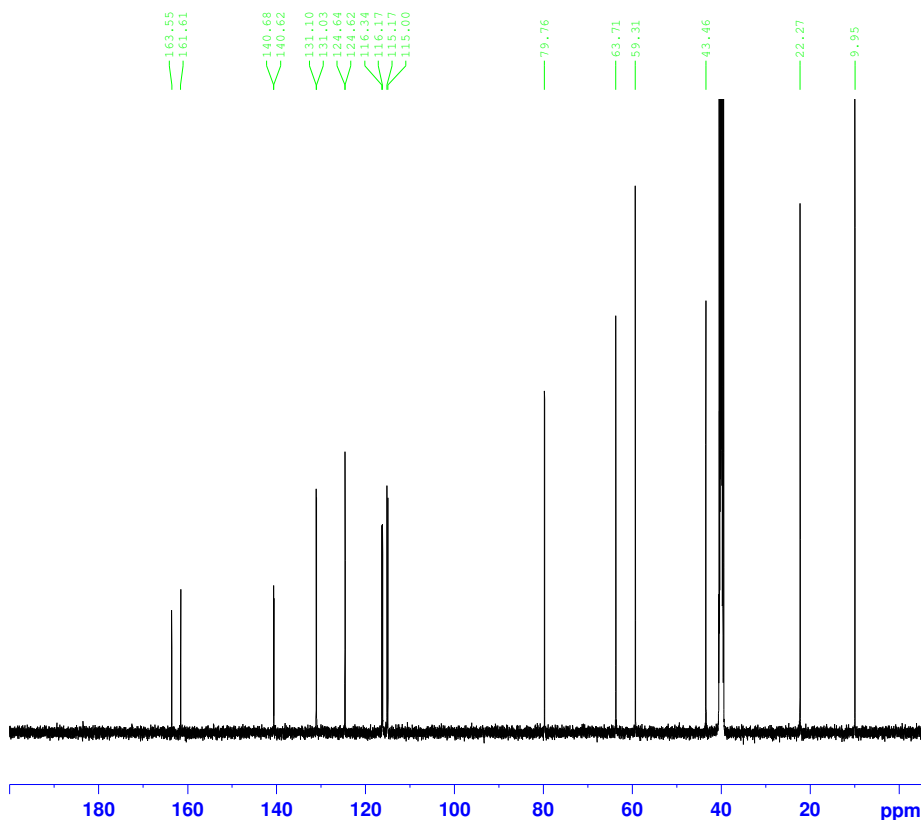

Current Data Parameters
NAME P-1292-15
EXPNO 2
PROCNO 1

F2 - Acquisition Parameters
Date_ 20151128
Time 23.43
INSTRUM spect
PROBHD 5 mm PABBO BB-
PULPROG zgpg30
TD 65536
SOLVENT DMSO
NS 3072
DS 4
SWH 29761.904 Hz
FIDRES 0.454131 Hz
AQ 1.1010048 sec
RG 2050
DW 16.800 usec
DE 6.50 usec
TE 296.0 K
D1 1.00000000 sec
D11 0.03000000 sec
TD0 1

===== CHANNEL f1 =====
SFO1 125.7703637 MHz
NUC1 13C
P1 9.00 usec
PLW1 122.00000000 W

===== CHANNEL f2 =====
SFO2 500.1320005 MHz
NUC2 1H
CPDPRG[2] waltz16
PCPD2 80.00 usec
PLW2 26.00000000 W
PLW12 0.32119001 W
PLW13 0.16186000 W

F2 - Processing parameters
SI 32768
SF 125.7577885 MHz
WDW EM
SSB 0
LB 1.00 Hz
GB 0
PC 1.40