




ANALYTICAL REPORT

5-MeO-N, N-DMT (C₁₃H₁₈N₂O)

2-(5-methoxy-1H-indol-3-yl)-N,N-dimethylethanamine

Remark – other NPS detected: **none**

Sample ID:	610-06A
Sample description:	powder - white
Sample type:	RM-reference material
Comments:	Sigma-Aldrich Lot#045K0956
Date of entry:	1/27/2006

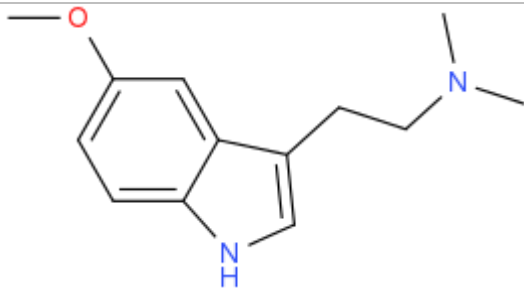Substance identified-structure ⁱ (base form)	
Systematic name	2-(5-methoxy-1H-indol-3-yl)-N,N-dimethylethanamine
Other names	5-methoxy-N,N-dimethyltryptamine
Formula (per base form)	C ₁₃ H ₁₈ N ₂ O
M _w (g/mol)	218,29
Salt form	
Smiles	COC=C1C=C2C(=CNC2=CC1)CCN(C)C
Compound Class	Indolalkylamines (fe tryptamines)
Other NPS detected	none
Add.info (purity..)	99%

This report has been produced with the financial support of the Prevention of and fight against crime Programme of the European Union (grant agreement number JUST/2013/ISEC/DRUGS/AG/6413). The contents of this report are the sole responsibility of the National Forensic Laboratory and can in no way be taken to reflect the views of the European Commission.

Supporting information

Analytical technique:	applied	remarks
GC-MS (EI ionization)	+	NFL GC-RT (min): 7,08 BP(1): 58; BP(2): 160,BP(3) :218,
FTIR-ATR	+	

GC-MS (Agilent): GC-method is RT locked to tetracosane (RT=9.53 min). Injection volume 1 ml and split mode (1:50) for GC-MS instruments and 1:5 for GC-MS-FTIR(condensed phase). Injector temperature: 280 °C. Column: HP1-MS (100% dimethylpolysiloxane), length 30 m, internal diameter 0.25 mm, film thickness 0.25 mm. Carrier gas He: flow-rate 1.2 ml/min. GC oven program: 170 °C for 1 min, followed by heating up to 293 °C at a rate of 18 °C/min, hold for 6.1 min, then heating at 50 °C/min up to 325 °C and finally 2.8 min isothermal. MSD source EI = 70 eV. GC-MS transfer line T= 235°C, source and quadropole temperatures 280°C and 180°C. m/z scan range: from 50 (40) to 550 amu.

FTIR-ATR (Perkin Elmer): scan range 4000-400 cm⁻¹; resolution 4cm⁻¹

FTIR (Spectra analyses-Danny): scan range 4000 to 700, resolution 4cm⁻¹

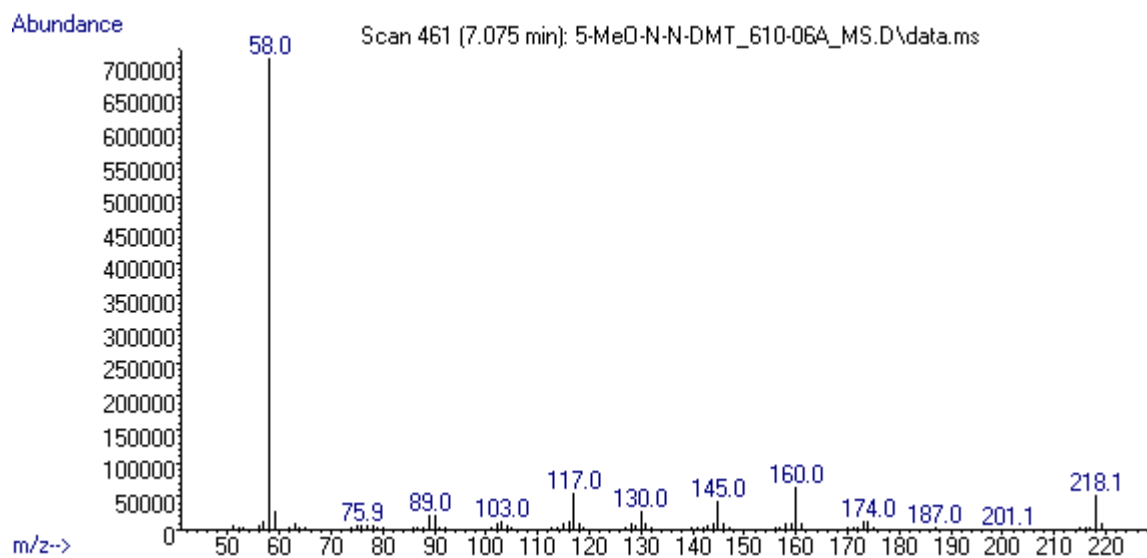


Figure 1: GC-MS spectrum



Figure 2: FTIR ATR

Created by OPSIN free tool: <http://opsin.ch.cam.ac.uk/> DOI: 10.1021/ci100384d