


## ANALYTICAL REPORT

AB-001 (C<sub>24</sub>H<sub>31</sub>NO)

**3-(adamantane-1-carbonyl)-1-pentyl-1H-indole**

Remark – other active cpd. detected: **none**

Sample ID:	1783-17
Sample description:	powder - white
Sample type:	RM-reference material
Comments <sup>1</sup> :	CHIRON Batch# 18 173; RESPONSE -purchasing
Date of entry:	3/31/2017

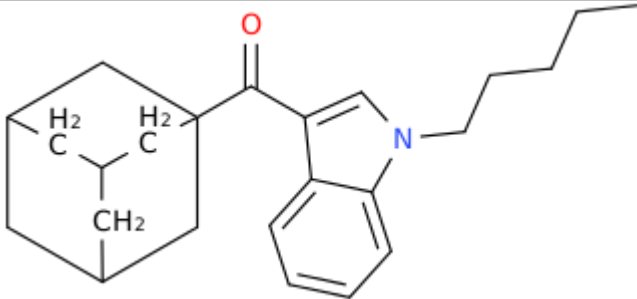Substance identified-structure <sup>2</sup> (base form)	
Systematic name:	3-(adamantane-1-carbonyl)-1-pentyl-1H-indole
Other names:	JWH-018 adamantoyl derivative; JWH-018 adamantyl analog; 1-pentyl-3-(1-adamantoyl)indole; 1-Adamantyl-(1-pentyl-1H-indol-3-yl)methanone
Formula (per base form)	C <sub>24</sub> H <sub>31</sub> NO
M <sub>w</sub> (g/mol)	349,52
Salt form:	base
StdInChIKey (for base form)	SHWDYCMMPWQM-UHFFFAOYSA-N
Other active cpd. detected	none
Add.info (purity..)	99%

<sup>1</sup> This report has been produced with the financial support of the Prevention of and fight against crime Programme of the European Union (grant agreement number JUST/2013/ISEC/DRUGS/AG/6413). The contents of this report are the sole responsibility of the National Forensic Laboratory and can in no way be taken to reflect the views of the European Commission.

<sup>2</sup> Created by OPSIN free tool: <http://opsin.ch.cam.ac.uk/> DOI: 10.1021/ci100384d



## Report updates

date	comments (explanation)

## Supporting information

Analytical technique:	applied	remarks
GC-MS (EI ionization)	+	NFL GC-RT (min): 14,14 BP(1): 214; BP(2): 144,BP(3) :215,
FTIR-ATR	+	direct measurement
GC-IR (condensed phase)	+	always as base form

**1. GC-MS** (Agilent): GC-method is RT locked to tetracosane (9.258 min). Injection volume 1 ml and split mode (1:50). Injector temperature: 280 °C. Chromatographic separation: on column HP1-MS (100% dimethylpolysiloxane), length 30 m, internal diameter 0.25 mm, film thickness 0.25 µm. Carrier gas He: flow-rate 1.2 ml/min. GC oven program: 170 °C for 1 min, followed by heating up to 190 °C at rate 8 °C/min, then heating up to 293 °C at a rate of 18 °C/min, hold for 7.1 min, then heating at 50 °C/min up to 325 °C and finally 6.1 min isothermal. MSD source EI = 70 eV. GC-MS transfer line T= 235°C, source and quadropole temperatures 280°C and 180°C, respectively. Scan range m/z scan range: from 50 (30 until 6 min.) to 550 (300 until 6 min) amu.

**2. FTIR-ATR** (Perkin Elmer): scan range 4000-400 cm<sup>-1</sup>; resolution 4cm<sup>-1</sup>

**3. GC- (MS)-IR** condensed phase (GC-MS (Agilent) & IR (Spectra analyses-Danny)

GC-method: Injection volume 1 ml and split mode (1:5). Injector temperature 280 °C. Chromatographic separation as above **(1)**. Split MS : IR = 1 : 9.

MSD source EI = 70 eV. GC-MS transfer line T= 235°C, source and quadropole temperatures 280°C and 180°C, respectively. Scan range m/z scan range: from 50 (30 until 6 min.) to 550 (300) amu.


IR (condensed (solid) phase): IR scan range 4000 to 650, resolution 4 cm<sup>-1</sup>.

4. HPLC-TOF for exact monoisotopic mass and empirical formula control - results are not shown in the report.


# ANALYTICAL RESULTS

MS (EI)

Abundance

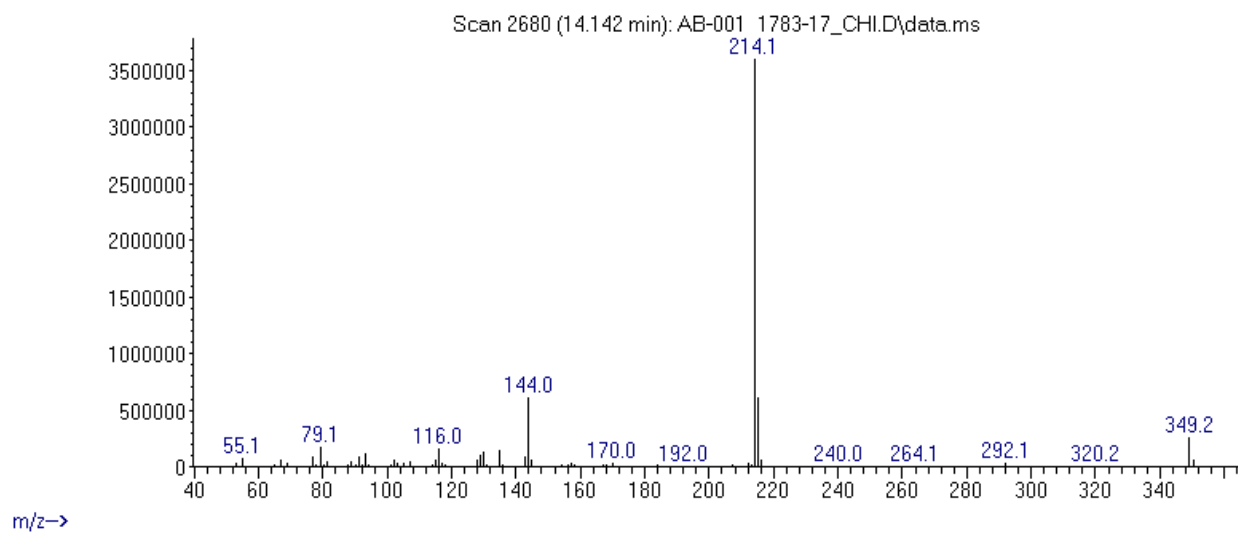

## FTIR-ATR - sample as received

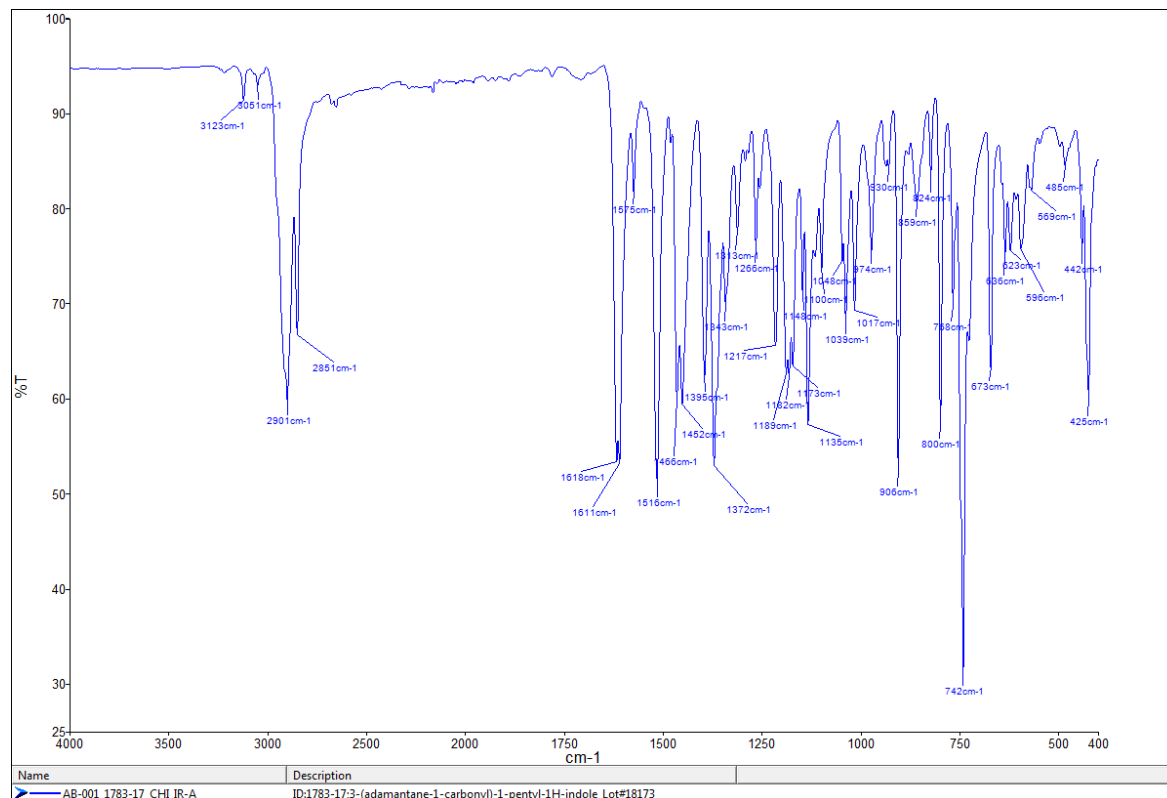

## IR (condensed phase – after chromatographic separation)

