


ANALYTICAL REPORT

DIPT (C₁₆H₂₄N₂)

N,N-Diisopropyltryptamine

Remark – other active cpd. detected: **none**

Sample ID:	1422-15
Sample description:	powder - white
Sample type:	RM-reference material
Comments ¹ :	Chiron AS Lot#16055; NFL- purchasing
Date of entry:	1/6/2016

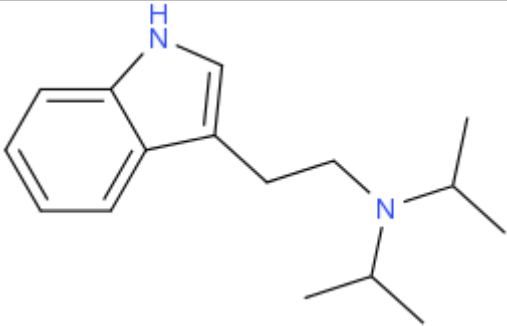Substance identified- structure ² (base form)	
Systematic name:	N,N-Diisopropyltryptamine
Other names:	
Formula (per base form)	C ₁₆ H ₂₄ N ₂
M _w (g/mol)	244.38
Salt form:	HCl
StdInChIKey	ZRVAAGAZUWXRIP-UHFFFAOYSA-N
Compound Class	Indolalkylamines (fe tryptamines)
Other active cpd. detected	none
Add.info (purity..)	99.20%

¹ This report has been produced with the financial support of the Prevention of and fight against crime Programme of the European Union (grant agreement number JUST/2013/ISEC/DRUGS/AG/6413). The contents of this report are the sole responsibility of the National Forensic Laboratory and can in no way be taken to reflect the views of the European Commission.

² Created by OPSIN free tool: <http://opsin.ch.cam.ac.uk/> DOI: 10.1021/ci100384d



Report updates

date	comments (explanation)

Supporting information

Analytical technique:	applied	remarks
GC-MS (EI ionization)	+	NFL GC-RT (min): 7.42 BP(1): 114; BP(2): 72, BP(3): 130,
FTIR-ATR	+	direct measurement
GC-IR (condensed phase)	+	

GC-MS (Agilent):

GC-method is RT locked to tetracosane (RT=9.53 min).

Injection volume 1 ml and split mode (1:50) .

Injector temperature: 280 °C.

Chromatographic separation

Column: HP1-MS (100% dimethylpolysiloxane), length 30 m, internal diameter 0.25 mm, film thickness 0.25 mm.

Carrier gas He: flow-rate 1.2 ml/min. GC oven program: 170 °C for 1 min, followed by heating up to 293 °C at a rate of 18 °C/min, hold for 6.1 min, then heating at 50 °C/min up to 325 °C and finally 2.8 min isothermal.

MSD source EI = 70 eV. GC-MS transfer line T= 235°C, source and quadrupole temperatures 280°C and 180°C, respectively. Scan range m/z scan range: from 50 (40) to 550 amu.

FTIR-ATR (Perkin Elmer): scan range 4000-400 cm⁻¹; resolution 4cm⁻¹

GC- (MS)-IR condensed phase (GC-MS (Agilent) & IR (Spectra analyses-Danny) IR scan range 4000 to 700, resolution 4cm⁻¹

GC-method:

Injection volume 1 ml and split mode (1:5) .

Injector temperature: 280 °C.

Chromatographic separation

Column: HP1-MS (100% dimethylpolysiloxane), length 30 m, internal diameter 0.25 mm, film thickness 0.25 mm.


Carrier gas He: flow-rate 1.2 ml/min. GC oven program: 170 °C for 1 min, followed by heating up to 293 °C at a rate of 18 °C/min, hold for 6.1 min, then heating at 50 °C/min up to 325 °C and finally 2.8 min isothermal.

Split MS : IR : (1:9)


MSD source EI = 70 eV. GC-MS transfer line T= 235°C, source and quadrupole temperatures 280°C and 180°C, respectively. Scan range m/z scan range: from 50 (40) to 550 amu.

IR (condensed phase): IR scan range 4000 to 700, resolution 4cm⁻¹

FIGURES OF SPECTRA

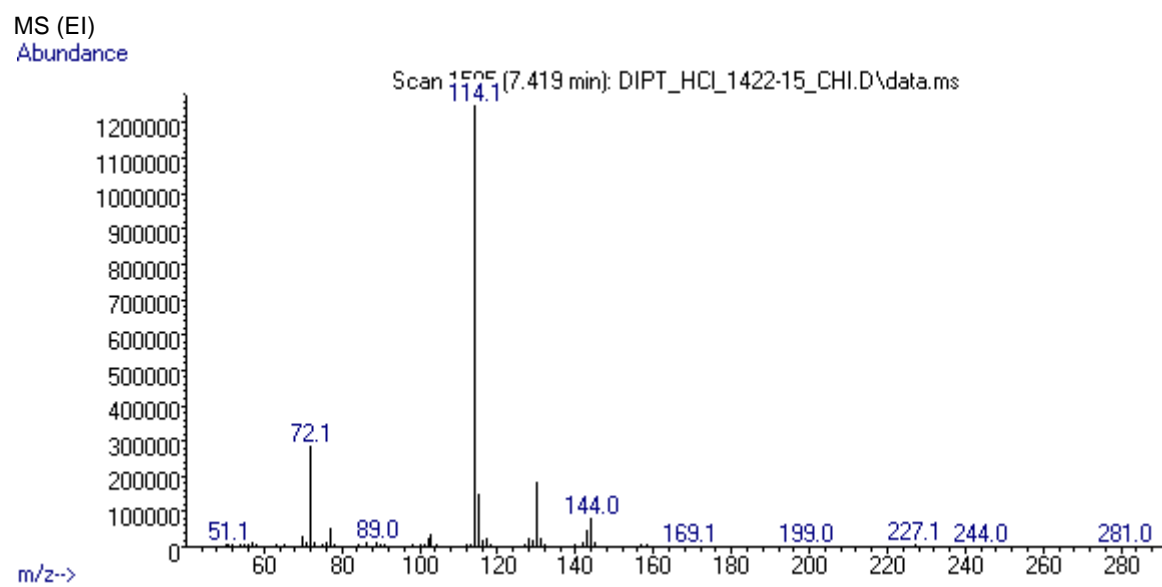

FTIR-ATR

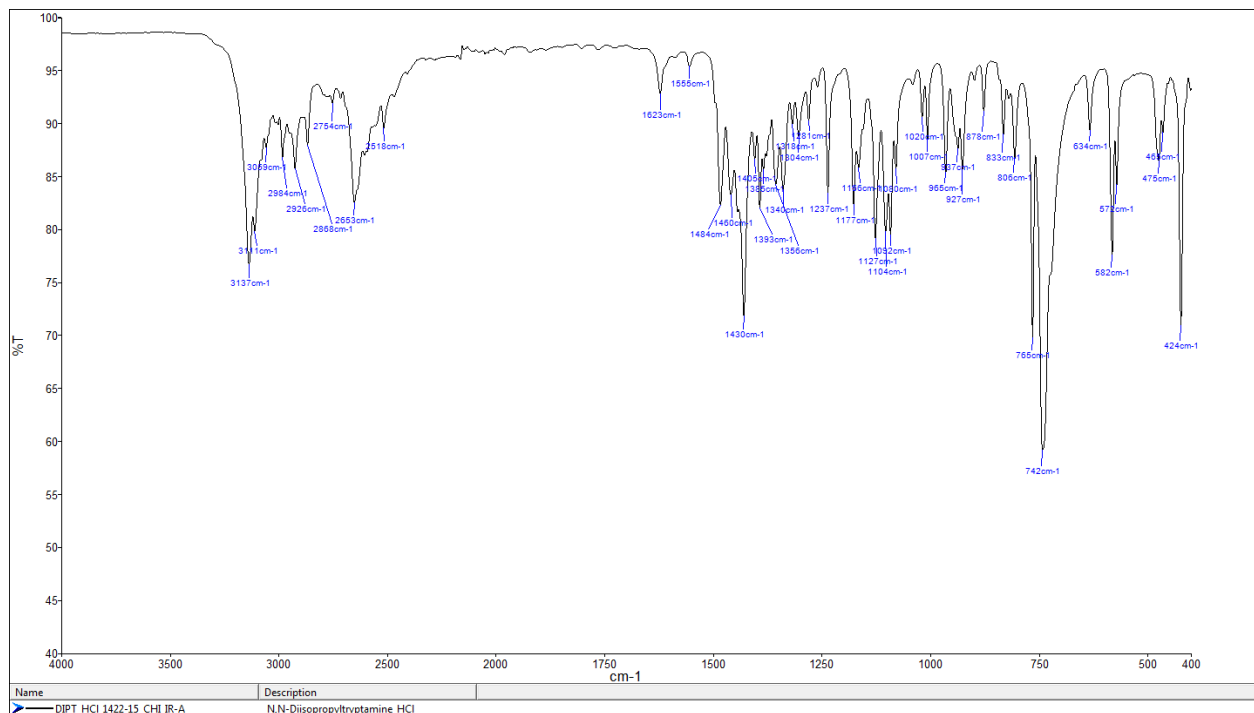

IR-Condensed phase

