


ANALYTICAL REPORT

Spirochlorphine (C₂₁H₂₄ClN₃O)

8-[1-(4-chlorophenyl)ethyl]-1-phenyl-1,3,8-triazaspiro[4.5]decan-4-one

Remark – other active cpd. detected: none

Sample ID:	3358-26
Sample description:	powder - white
Sample type:	RM-reference material
Comments:	CAY Lot#0809267-6,
Date of entry (DD/MM/YYYY):	27/08/2026
web link	http://www.policija.si/apps/nfl_response_web/seznam.php

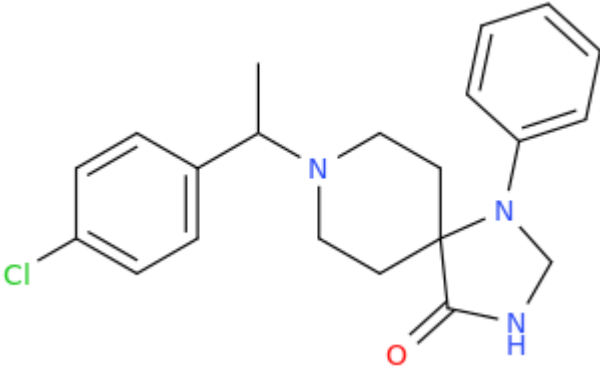Substance identified- structure ¹ (base form)	
Systematic name:	8-[1-(4-chlorophenyl)ethyl]-1-phenyl-1,3,8-triazaspiro[4.5]decan-4-one
Other names:	R-6890; 8-(p-chloro-α-methylbenzyl)-1-phenyl-1,3,8-triazaspiro[4.5]decan-4-one; 1-phenyl-4-oxo-8-[1-(4-chloro-phenyl)-ethyl]-1,3,8-triaza-spiro(4.5)decane; 8-[1-(4-chlorophenyl)-ethyl]-4-oxo-1-phenyl-1,3,8-triazaspiro[4.5]decane; 1-phenyl-4-oxo-8-[1-(
Formula (per base form)	C ₂₁ H ₂₄ ClN ₃ O
M _w (g/mol)	369,89
Salt form:	base
StdInChIKey (per base form)	KFEYPBZJPJRF-XUHFFFAOYSA-N
Other active cpd. detected	none
Add.info (purity..)	≥98%

¹ Created by OPSIN free tool: <http://opsin.ch.cam.ac.uk/> DOI: 10.1021/ci100384d

Report updates

date	comments (explanation)

Supporting information

Analytical technique:	applied	remarks
GC-MS (EI ionization)	+	NFL GC-RT (min): 15,53 BP(1): 175; BP(2): 139,BP(3) :56,
FTIR-ATR	+	direct measurement
GC-IR (solid phase)	+	always as base form
HPLC-TOF	+	exact mass theoretical: 369,1608 / measured Δ ppm: -2,48

1. GC-MS (Agilent): GC-method is RT locked to tetracosane (9.258 min). Injection volume 1 μ l and split mode (1:50). Injector temperature: 280 °C. Chromatographic separation: on column HP1-MS (100% dimethylpolysiloxane), length 30 m, internal diameter 0.25 mm, film thickness 0.25 μ m. Carrier gas He: flow-rate 1.2 ml/min. GC oven program: 170 °C for 1 min, followed by heating up to 190 °C at rate 8 °C/min, then heating up to 293 °C at a rate of 18 °C/min, hold for 6.1 min, then heating at 50 °C/min up to 325 °C and finally 9.1 min isothermal. MSD source EI = 70 eV. GC-MS transfer line T= 235°C, source and quadropole temperatures 280°C and 180°C, respectively. Scan range m/z scan range: from 50 (30 until 6 min.) to 550 (300 until 6 min) amu.

2. FTIR-ATR (Perkin Elmer): scan range 4000-400 cm^{-1} ; resolution 4 cm^{-1}

3. GC- (MS)-IR solid phase (GC-MS (Agilent) & IR (Spectra analyses-Danny)

GC-method: Injection volume 1 μ l and split mode (1:5). Injector temperature 280 °C. Chromatographic separation as above **(1)**. Split MS: IR = 1 : 9.

MSD source EI = 70 eV. GC-MS transfer line T= 235°C, source and quadropole temperatures 280°C and 180°C, respectively. Scan range m/z scan range: from 50 (30 until 6 min.) to 550 (300) amu.


IR (condensed (solid) phase): IR scan range 4000 to 650, resolution 4 cm^{-1} .

4. HPLC-TOF (Agilent): 6230B TOF with Agilent 1260 Infinity HPLC with binary pump, column: Zorbax Eclipse XDB-C18, 50 x 4.6 mm, 1.8 micron. Mobile phases (A) 0.1% formic acid and 1mM ammonium formate in water; (B) 0.1% formic acid in methanol (B). Gradient: starting at 5% B, changing to 40% B over 4 min, then to 70% over 2 min and in 5 min to 100%, hold 1 min and back to 5%, equilibration for 1.7 min. The flow rate: 1.0 ml/min; Injection volume 1 μ l. MS parameters: 2GHz, Extended Dynamic range mode to a maximum of 1700 amu, acquisition rate 1.30 spectra/sec. Sample ionisation: by Agilent Jet Stream technology (Dual AJS ESI). Ion source: positive ion scan mode with mass scanning from 82 to 1000 amu. Other TOF parameters: drying gas (N₂) and sheath temperature 325 °C; drying gas flow rate 6 l/min; sheath gas flow rate 8 l/min; nebulizer 25 psig; Vcap. 4000 V; nozzle 2000 V; skimmer 65 V; fragmentor 175 V and Octopole RF 750 V.

ANALYTICAL RESULTS

MS (EI)

Abundance

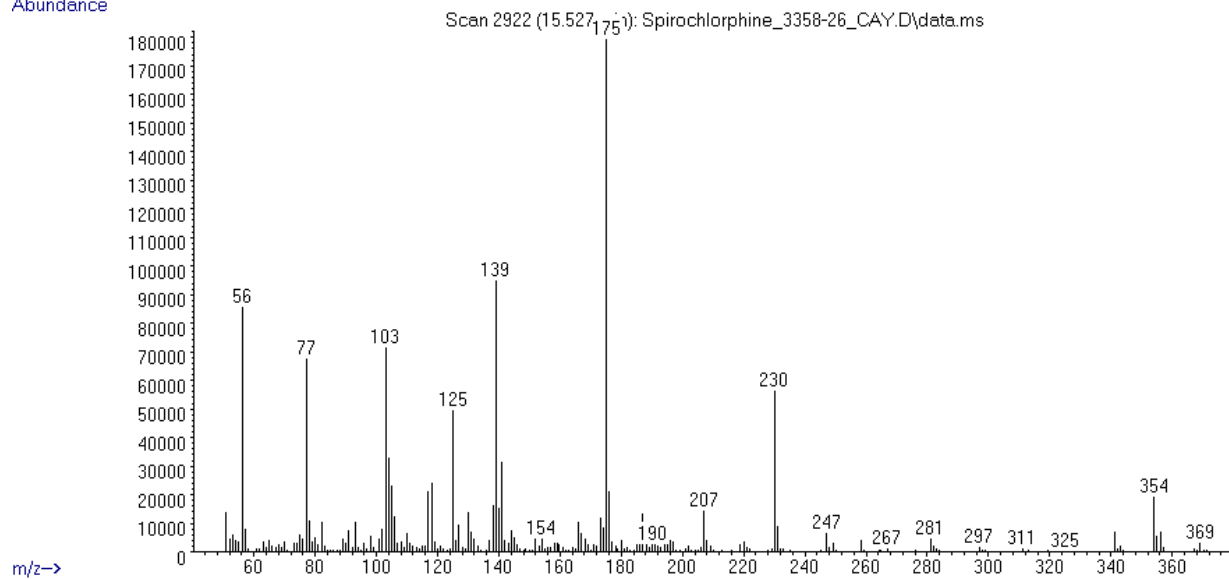

Name
 > Spirochlorphine 3358-26, CAV, IR-A

Description
 ID: 3358-26; 8-[1-[4-chlorophenyl]ethyl]-1-phenyl-1,3,8-triazaspiro[4.5]decan-4-one, Lot#0809267-6

Peak Data (Wavenumbers in cm⁻¹):

Wavenumber (cm ⁻¹)
3067.99
2843.45
1704.07
1597.63
1573.70
1501.68
1472.34
1387.89
1369.93
1348.64
1301.67
1187.80
1157.17
1119.95
1111.56
1093.46
1033.58
1011.56
940.79
933.35
933.68
878.70
836.34
825.11
807.82
787.70
740.72
691.66
648.45
611.48
544.80
517.19
482.82
442.51

NOTE: This is condensed phase IR (per base form of substance)
Instrument (DiscovIR-GC)

Peak Wavenumber (cm ⁻¹)
3212.36
3103.21
3062.30
2930.35
2971.22
2829.45
1708.38
1601.06
1571.16
1503.16
1470.51
1392.55
1348.71
1322.65
1196.08
1119.38
1092.01
1014.13
923.56
833.31
748.49
742.70
714.68
694.96