




ANALYTICAL REPORT¹

TH-PVP (C₁₉H₂₇NO)

2-(pyrrolidin-1-yl)-1-(5,6,7,8-tetrahydronaphthalen-2-yl)pentan-1-one

Remark – other NPS detected: **none**

Sample ID:	1278-15
Sample description:	crystalinic - beige (gently)
Sample type:	test purchase /RESPONSE -purchasing
Date of sample receipt (M/D/Y):	9/18/2015
Date of entry (M/D/Y) into NFL database:	10/4/2015
Report updates (if any) will be published here:	http://www.policija.si/apps/nfl_response_web/seznam.php

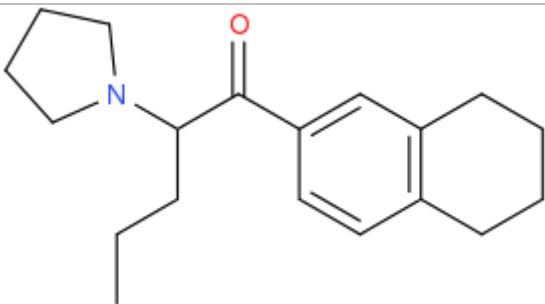Substance identified - structure ² (base form)	
Systematic name	2-(pyrrolidin-1-yl)-1-(5,6,7,8-tetrahydronaphthalen-2-yl)pentan-1-one
Other names	
Formula (per base form)	C ₁₉ H ₂₇ NO
M _w (g/mol)	285,43
Salt form	HCl
StdInChIKey	MMIKQWIZKBYLKZ-UHFFFAOYSA-N, MMIKQWIZKBYLKZ-UHFFFAOYSA-N
Compound Class	Cathinones
Other NPS detected	none
Add.info (purity..)	pure by GC-MS, HPLC-TOF, pure by NMR

¹ This report has been produced with the financial support of the Prevention of and fight against crime Programme of the European Union (grant agreement number JUST/2013/ISEC/DRUGS/AG/6413). The contents of this report are the sole responsibility of the National Forensic Laboratory and can in no way be taken to reflect the views of the European Commission.

² Created by OPSIN free tool: <http://opsin.ch.cam.ac.uk/> DOI: 10.1021/ci100384d

Report updates

date	comments (explanation)
08/12/2016	substance class changed

Instrumental methods (if applied) in NFL

1. GC-MS (Agilent): GC-method is RT locked to tetracosane (RT=9.53 min). Injection volume 1 ml and split mode (1:50) . Injector temperature: 280 °C. Chromatographic separation: on column HP1-MS (100% dimethylpolysiloxane), length 30 m, internal diameter 0.25 mm, film thickness 0.25 mm. Carrier gas He: flow-rate 1.2 ml/min. GC oven program: 170 °C for 1 min, followed by heating up to 293 °C at a rate of 18 °C/min, hold for 6.1 min, then heating at 50 °C/min up to 325 °C and finally 2.8 min isothermal. MSD source EI = 70 eV. GC-MS transfer line T= 235°C, source and quadropole temperatures 280°C and 180°C, respectively. Scan range m/z scan range: from 50 (40) to 550 amu.

2. HPLC-TOF (Agilent): 6230B TOF with Agilent 1260 Infinity HPLC with binary pump, column: Zorbax Eclipse XDB-C18, 50 x 4.6 mm, 1.8 micron. Mobile phases (A) 0.1% formic acid and 1mM ammonium formate in water; (B) 0.1% formic acid in methanol (B). Gradient: starting at 5% B, changing to 40% B over 4 min, then to 70% over 2 min and in 5 min to 100%, hold 1 min and back to 5%, equilibration for 1.7 min. The flow rate: 1.0 ml/min; Injection volume 1 µl. MS parameters: 2GHz, Extended Dynamic range mode to a maximum of 1700 amu, acquisition rate 1.30 spectra/sec. Sample ionisation: by Agilent Jet Stream technology (Dual AJS ESI). Ion source: positive ion scan mode with mass scanning from 82 to 1000 amu. Other TOF parameters: drying gas (N₂) and sheath temperature 325 °C; drying gas flow rate 6 l/min; sheath gas flow rate 8 l/min; nebulizer 25 psig; Vcap. 4000 V; nozzle 2000 V; skimmer 65 V; fragmentor 175 V and Octopole RF 750 V.

3. FTIR-ATR (Perkin Elmer): scan range 4000-400 cm⁻¹; resolution 4cm⁻¹

4. GC- (MS)-IR condensed phase (GC-MS (Agilent) & IR (Spectra analyses-Danny)

GC-method: Injection volume 1 ml and split mode (1:5). Injector temperature 280 °C. Chromatographic separation as above (**1**). Split MS : IR = 1:9.

MSD source EI = 70 eV. GC-MS transfer line T= 235°C, source and quadropole temperatures 280°C and 180°C, respectively. Scan range m/z scan range: from 50 (40) to 550 amu.

IR (condensed phase): IR scan range 4000 to 650, resolution 4 cm⁻¹.

5. IC (anions) (Thermo Scientific, Dionex ICS 2100), Column: IonPac AS19, 2 x 250mm; Eluent: 10mM from 0 to 10 min, 10-58 mM from 10 to 40min; Flow rate: 0.25 ml/min; Temperature: 30°C; Suppressor: AERS 500 2mm, suppressor current 13mA; Inj. Volume: 25 µl

Supporting information


Solubility in	result/remark
CH ₂ Cl ₂	partially
MeOH	soluble
H ₂ O	

Analytical technique:	applied	remarks
GC-MS (EI ionization)	+	NFL GC-RT (min): 8,79 BP(1): 126; BP(2): 127,BP(3) :91,
HPLC-TOF	+	Exact mass (theoretical): 285,2093; measured value Δppm:-0,53; formula:C19H27NO
FTIR-ATR	+	direct measurement (sample as received)
FTIR (condensed phase) always as base form	+	
IC (anions)	+	
NMR (in FKKT)	+	
validation		
other		


ANALYTICAL RESULTS

MS (EI)


Abundance

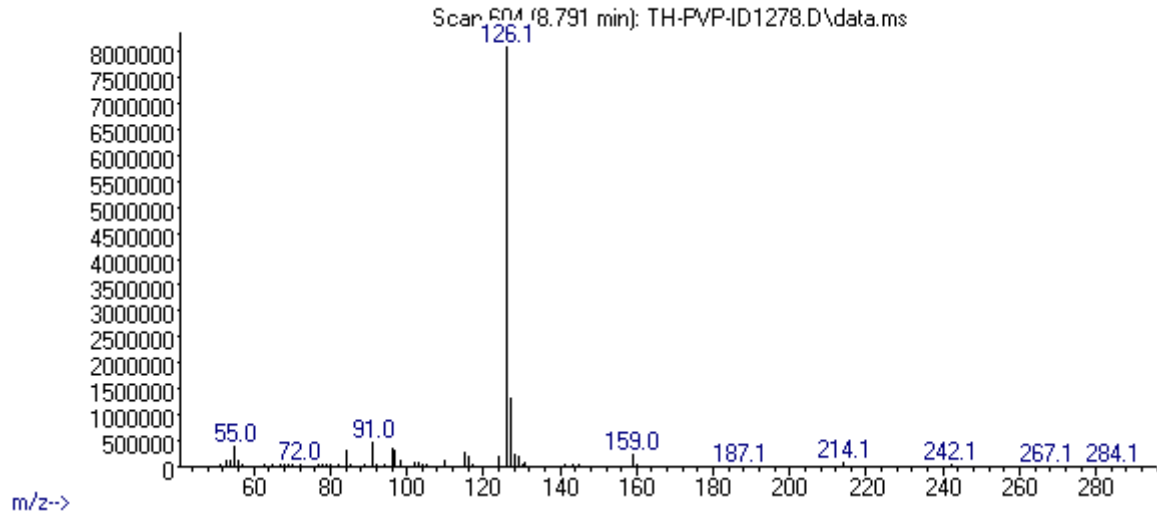

FTIR-ATR - direct measurement (sample as received)

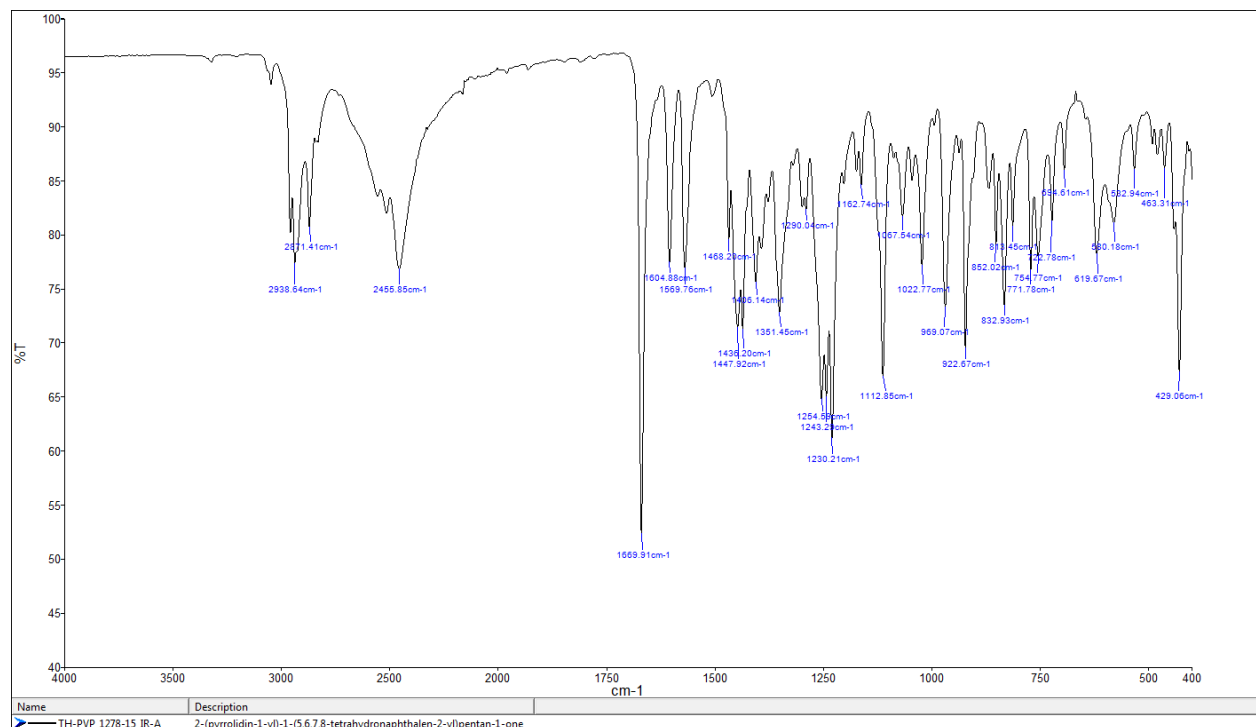

IR (condensed phase – after chromatographic separation)

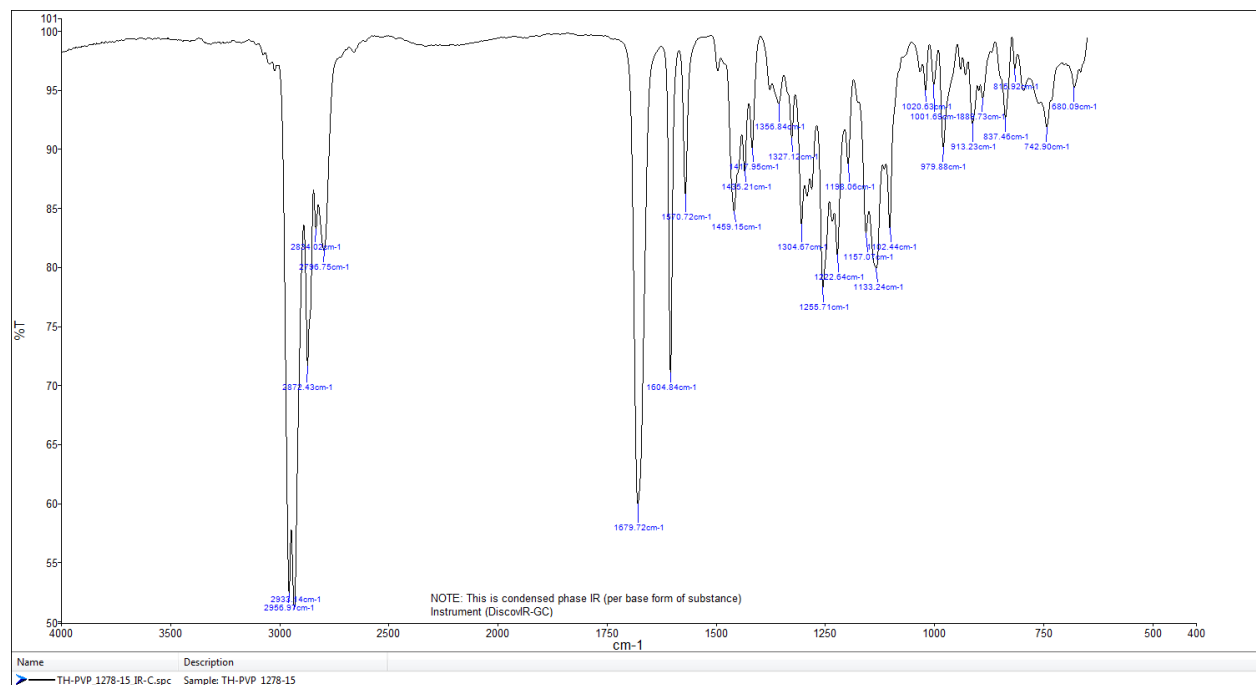

TOF REPORT


Data File	TH-PVP_1278-15_TOF.d	Sample Name	TH-PVP
Sample Type	Sample	Position	P1-E3
Instrument Name	6230B TOF LC-MS	User Name	TG
Acq Method	droge general-13-5-2015-XDB-C18-ESI-poz.m	Acquired Time	9/21/2015 9:43:26 AM
IRM Calibration Status	Success	DA Method	Droge_Default.m
Comment	extract in MeOH		

Compound Table


Label	Compound Name	Obs. RT	Obs. Mass
Cpd 1: TH-PVP	TH-PVP	7.277	285.2094

Name	Obs. m/z	Obs. RT	Obs. Mass	DB RT	DB Formula	DB Mass	DB Mass Error (ppm)
TH-PVP	286.2167	7.277	285.2094	7.3	C19 H27 N O	285.2093	-0.53


Compound Chromatograms

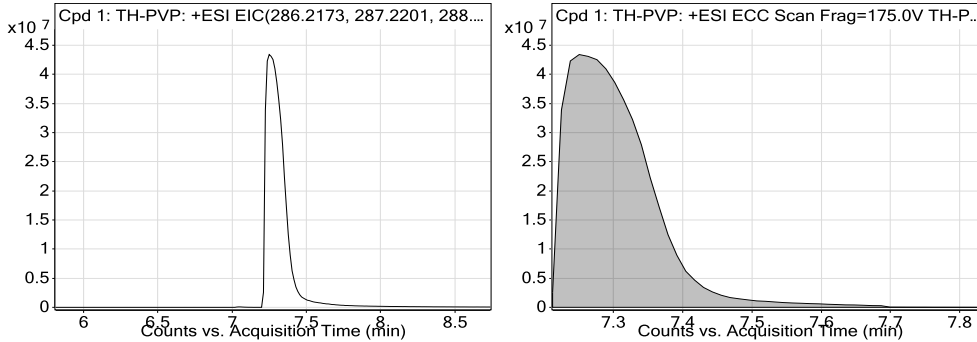

MFE MS Zoomed Spectrum

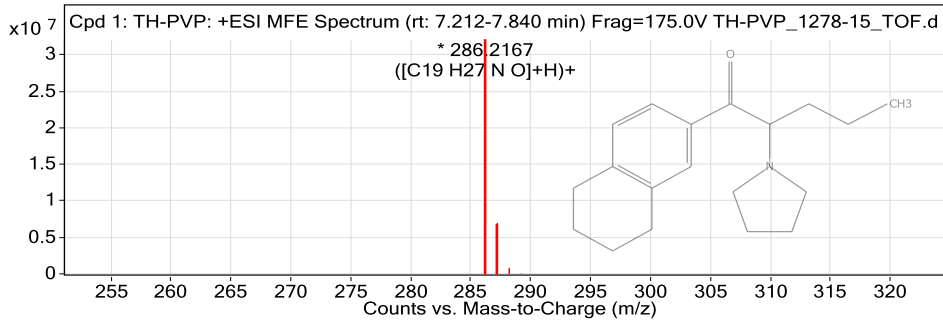

MS Zoomed Spectrum

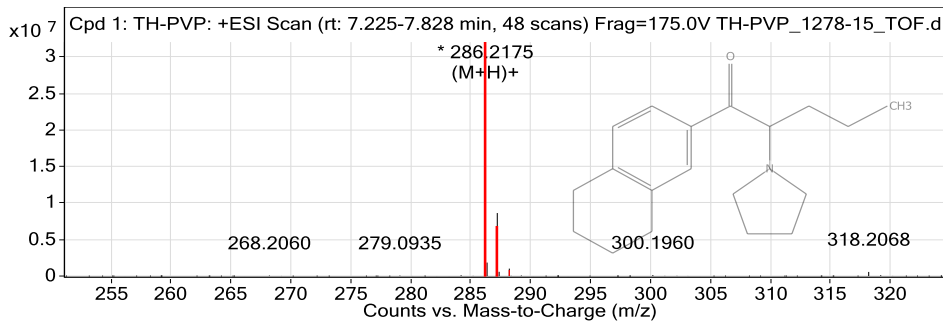

MS Spectrum Peak List


Obs. m/z	Charge	Abund	Formula	Ion/Isotope
286.2167	1	32040584	C19 H27 N O	(M+H)+
287.22	1	6946613.01	C19 H27 N O	(M+H)+
288.2236	1	681537.96	C19 H27 N O	(M+H)+
289.2265	1	49168.92	C19 H27 N O	(M+H)+
290.228	1	1332.51	C19 H27 N O	(M+H)+

--- End Of Report ---

Peak Integration Report

Sample Name:	TH-PVP_1278-15_IC	Inj. Vol.:	25,00
Injection Type:	Unknown	Dilution Factor:	1,0000
Program:	ANIONI	Operator:	kemija
Inj. Date / Time:	27-okt-2015 / 12:16	Run Time:	42,00

No.	Time min	Peak Name	Peak Type	Area $\mu\text{S} \cdot \text{min}$	Height μS	Amount mg/L
1,00	9,48	Chloride	BMB	3,23	13,59	n.a.
TOTAL:				3,23	13,59	0,00

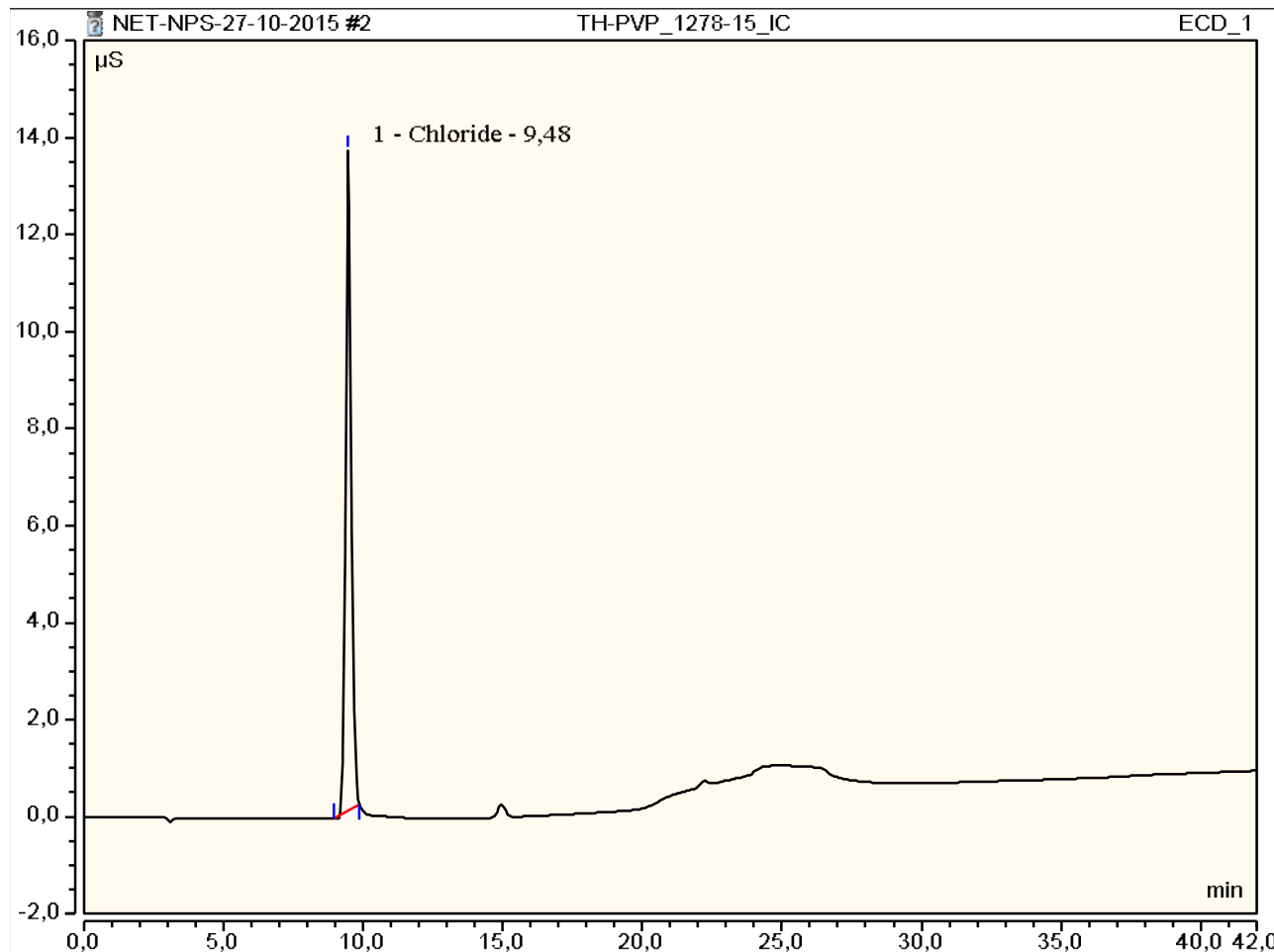



REPORT

Sample ID:	1278-15
Our notebook code:	P-1278-15
NMR sample preparation:	15 mg dissolved in 0.7 mL DMSO- d_6
NMR experiments:	^1H , ^{13}C .
Proposed structure:	
Chemical name:	1-(1-oxo-1-(5,6,7,8-tetrahydronaphthalen-2-yl)pentan-2-yl)pyrrolidin-1-ium
Comments:	- Structure elucidation based on 1D NMR spectra - Compound is pure by NMR.
Supporting information:	Copies of ^1H and ^{13}C NMR spectra
Author:	Prof. Dr. Janez Košmrlj, Doc. Dr. Krištof Kranjc
Date of report:	December 8, 2015

P-1278-15



```


Current Data Parameters
NAME          P-1278-15
EXPNO         1
PROCNO        1

F2 - Acquisition Parameters
Date_         20151129
Time          4.53
INSTRUM       spect
PROBHD        5 mm PABBO BB-
PULPROG       zg30
TD            65536
SOLVENT       DMSO
NS            16
DS            2
SWH           10000.000 Hz
FIDRES        0.152588 Hz
AQ            3.2768500 sec
RG            114
DW            50.000 usec
DE            6.50 usec
TE            296.0 K
D1            1.00000000 sec
TD0           1

===== CHANNEL f1 =====
SFO1          500.1330885 MHz
NUC1           1H
P1             8.90 usec
PLW1          26.00000000 W

F2 - Processing parameters
SI            65536
SF            500.1300024 MHz
WDW           EM
SSB           0
LB            0.30 Hz
GB            0
PC            1.00

```

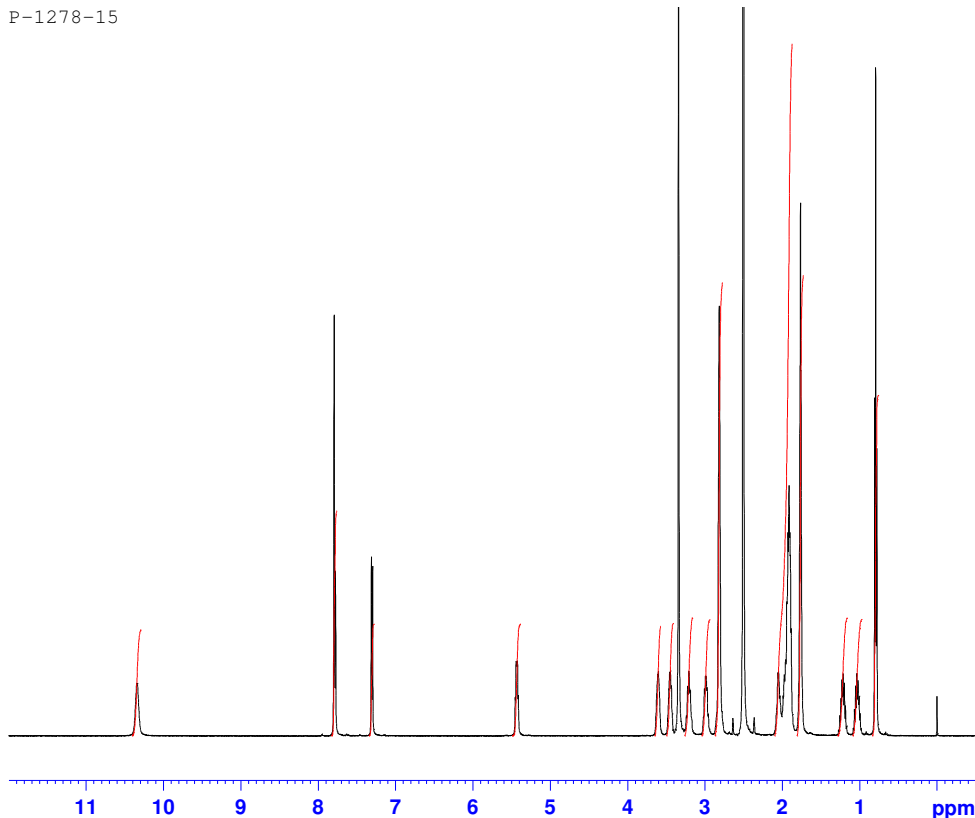

P-1278-15



```

Current Data Parameters
NAME          P-1278-15
EXPNO         2
PROCNO        1

F2 - Acquisition Parameters
Date_         20151129
Time          7.22
INSTRUM       spect
PROBHD        5 mm PABBO BB-
PULPROG       zgpg30
TD            65536
SOLVENT       DMSO
NS            4096
DS            4
SWH           29761.904 Hz
FIDRES        0.454131 Hz
AQ            1.1010048 sec
RG            2050
DW            16.800 usec
DE            6.50 usec
TE            296.0 K
D1            1.00000000 sec
D11           0.03000000 sec
TD0           1

===== CHANNEL f1 =====
SFO1          125.7703637 MHz
NUC1           13C
P1             9.00 usec
PLW1          122.00000000 W

===== CHANNEL f2 =====
SFO2          500.1320005 MHz
NUC2           1H
CPDPRG[2]     waltz16
PCPD2         80.00 usec
PLW2          26.00000000 W
PLW12         0.32179001 W
PLW13         0.16186000 W

F2 - Processing parameters
SI            32768
SF            125.7577885 MHz
WDW           EM
SSB           0
LB            1.00 Hz
GB            0
PC            1.40

```

